

2024

Aseptico

GLOBAL DENTAL CATALOG

IMPLANT MOTORS

NEW MOTOR: 1070 Series (pg. 4)

DENTAL INNOVATIONS

Monitor Osseointegration: Penguin II (pg. 22)

PORTABLE DENTISTRY

Aseptico's Most Advanced Dental Unit:
GO Ultra Portable Dental System (pg. 34)



SONYA DUNBAR MPH, RDH
Mobile Dental Xpress
Jacksonville, Florida

EXCLUSIVE COUPON INSIDE

Welcome to the 2024 Aseptico Global Dental Catalog

This year, we are proud to introduce our brand-new AEU-1070 implant motor. Designed for implant focused dentists, it is the ultimate tool for simplifying workflow. It makes a fine addition to our motor line as an efficient “implants only” system and continues to prove our commitment to providing our customers with the most innovative products. If you place or restore implants, the AEU-1070 implant motor will make an exceptional addition to your practice. See it featured on page 4.

We are pleased to report that the GO Ultra Portable Dental System continues to excel in a variety of settings — including urban, rural, and humanitarian clinics. Users rave about the lightweight system and how quiet the unit is. Check out the link on page 34 about how Sonya Dunbar utilizes her GO Ultra Portable dental system and how it revolutionized her day-to-day work.

The products of which we are most proud, Aseptico dental motors, have earned numerous industry accolades and a worldwide reputation for legendary performance and reliability. In these pages you will find our complete line, including the award-winning 7000 Series which has been named both a Top 25 Endodontic Product and a Top 25 Dental Implant Product by Dentistry Today magazine for 9 consecutive years.

We take pride that almost all Aseptico products are made in the USA and invite you to compare our products to any other for quality, performance, innovative features, reliability and value.



Glenn Kazen, President



20% OFF ENTIRE ORDER

Use Code **Catalog24** to redeem this Valuable Offer when you call or purchase online!

Coupon valid January 1, 2024 – April 30, 2024. Some exclusions apply*

*Take 20% off entire order between January 1, 2024 and April 30, 2024. Some exclusions apply: the AMC-25 Self-Contained Mobile Operator, Transport III Portable Dental System (AEU-525/AEU-525S), NOMAD Portable X-ray Systems, HandiDam® Pre-framed Dental Dams, Spotter Smart Implant Detector, NSK handpieces, W&H handpieces and W&H handpiece lubricant are excluded. Valid for one use only. This offer cannot be combined with any other offer. Valid only on new orders, NOT repairs. Some other exclusions might apply. Please call our friendly and professional team with any questions at 425-408-9347.

The Aseptico Global Dental Catalog is published annually by Aseptico Inc., 8333 216th Street SE, Woodinville, WA 98072. Aseptico is a global leader in dental equipment design and engineering. The company has endeavored to represent all current product information, specifications and pricing accurately. Product information may change. If you would like to communicate with Aseptico, the company welcomes your interest and invites you to contact us at 425-408-9347 or via asepticocom.com. © Aseptico Inc., 2023 all rights reserved. *Aseptico® is a registered Trademark of Aseptico Inc.

Aseptico

GLOBAL DENTAL CATALOG PROVEN PRODUCTS



Dental Motors & Handpieces WE POWER DENTISTRY

Oral Surgery and Implant Motors	3-8
NEW 1070 Series Implant Motor	4
Electrical Operative Motor	9
Endodontic Motors	10-11
Handpieces	12-18
Tubing, Accessories & Replacement Parts	19-20



Innovations & Utility PROVEN PRODUCTS

Penguin II Implant Stability Monitor	22
Spotter Smart Implant Detector	27
HandiDam® Pre-Framed Dental Dam	28
Ortho Technique Chain & Grid	28
Dental Dam Clamps	28
Dental Dam Punches and Forceps	28
Dental Delivery Systems	29-30
Powered Trolley Carts	31
NEW Deluxe Operatory Cart	32



Mobile & Portable Dentistry MOVING TO MEET NEED

Portable Delivery Systems	34-39
GO Ultra Portable Delivery System	30-31
Mobile Dental Carts	40-41
Portable Chairs & Stools	42-43
NEW NOMAD Pro Portable X-ray	44
Portable LED Exam Light	45
Portable Tray Stand	45
Lab Motors	46
Portable Saliva Ejector	47

Customer Service & Ordering VALUE & EXCELLENCE

Ordering Information, Policies and Contact Information	48
Repair and Product Support	48

Call 425.408.9347 or Shop Online at www.aseptico.com

Aseptico® DENTAL MOTORS AND HANDPIECES



2 YEAR WARRANTY
on Aseptico® Dental Motors Made in the USA*

*1 Year warranty on handpieces. Read our warranty policy on page 48



Find The Right Dental Motor For You

All Aseptico® Dental Motors and Micromotors are made in the USA. Each surgical system features autoclavable tubing which saves money and is better for the planet. With a variety of different motor options available, Aseptico makes it easy to find the perfect pairing for your practice. See the table below for features.

Motor Features	AEU-7000	AEU-6000	AEU-1070	AEU-27A
Oral Surgery	✓	✓		
Restorative	✓	✓		
Implantology	Up to 80 Ncm	Up to 50 Ncm	Up to 70 Ncm	
Rotary Endodontics	✓	✓		✓
Auto-Stop Reverse for Endo	✓	✓		✓
Dynamometer	✓			
Upgradeable Software	✓	✓		✓
Programmable Presets	Up to 6	Up to 5	Up to 8	File Library*
Multi-Function Foot Control	✓	✓	✓	
Autoclavable Tubing	✓	✓	✓	

*Dentsply Sirona File Library Included

"We were in need of a reliable implant/oral surgery motor for our busy practice that includes an Oral Surgeon. With such a heavy motor workload daily, it was imperative to pick a system with enough cutting power for both procedures daily. The 7000 Series is one piece of equipment that I would not want to practice without."



Erich Jurges, DDS
Sage Brush Dental

"Aseptico® surgical and implant motors have been a staple in my Oral and Maxillofacial Surgery practice. The quality, durability, ease of use and customer service have been unmatched. I wholeheartedly recommend them for all your surgical motor and implant handpiece needs."



Joyce Jeffries, DDS
Cottonwood Oral & Maxillofacial Surgery

CALL TODAY! 425.408.9347 or shop online at aseptico.com

1070 Series Our Newest Implant Motor

Streamline Your Implant Procedures

AEU-1070-70V



"It's amazing that it is such a workhorse! I've placed over 40,000 implants and have had practically no issues with these machines! I also very much like the built in irrigation and the convenient multi-function foot pedal!"



Todd E. Shatkin, DDS, DICOI, DIAMDI, FICD, FADI, FIADI
Shatkin First

Streamlined workflow

The 1070 Series Implant motor has the ability to place implants up to 70 Ncm with a 20:1 reduction handpiece and features a 1:1 setting for surgical and low-speed procedures.

User-friendly interface

The touchscreen display and intuitive user interface make it simple to adjust motor settings based on your preferences.

Hands-free control

The multifunction foot control enables easy adjustment of speed, direction, irrigation flow, torque, and presets, enhancing efficiency and infection control.

Precise and powerful

The brushless 40,000 RPM motor delivers the power needed for accurate and efficient implant procedures.

Versatility

With the 1:1 setting (and available 1:1 straight handpiece attachment) the motor is also suitable for surgical and low-speed procedures as well.

Customization options

The eight programmable presets allow you to customize the motor settings based on your preferred techniques.

Enhanced accuracy

The integrated handpiece calibration feature of the 1070 Series Implant Motor ensures speed and torque accuracy at the bur, providing consistent and reliable performance.

Compact and durable

The motor's compact footprint and water-resistant chassis design make it suitable for various clinical environments.

Sustainability-focused

The 1070 Series Implant Motor uses autoclavable irrigation tubing. That helps in reducing waste and increasing saving costs compared to single-use, disposable tubing.

Watch a Video
Demonstration



Blog:
FAQ on the new 1070
Series Implant Motor



1070 Motor with On-Off Foot Switch	AEU-1070	\$3,330
1070 Motor with Multi-Function Foot Control	AEU-1070-70V	\$3,550

Options

Micromotor w/ 10ft Cable	AE-230A-40EXT	\$1,200
Multi-Function Foot Control w/ 10ft Cable	AE-70V2-10	\$660
On/Off Foot Control w/ 10ft Cable	AE-7PM-10	\$210

Replacement Parts

Micromotor w/ 6ft Cable	AE-230A-40	\$1,100
Multi-Function Foot Control	AE-70V2	\$610
On/Off Foot Control	AE-7PM	\$190
Irrigation Bag Hanger Rod	461541	\$48
Power Cord	840107	\$23
Handpiece Cradle	461561-01	\$11
Irrigation Tubing Set	AE-24	\$80

PACKAGE DISCOUNTS AVAILABLE.
CALL TODAY! 425.408.9347

EXCLUSIVE 1070 PROMO ON BACK OF CATALOG

Mont Blanc AHP-85MB-X shown (not included). See details and other compatible handpieces pages 12-18.

7000 Series Motor All-in-One Surgical Platform

Oral Surgery, Implant Placement,
Endodontics & Restorative Dentistry

AEU-7000



Mont Blanc AHP-85MBFO-CX shown (not included).
See details and other compatible handpieces pages 12-18.

The Aseptico 7000 Series Dental Motor

is one of the most powerful and easy-to-use implant/surgical motors available today. This all-in-one platform provides specialist-level performance for restorative, oral surgery, implant, and rotary endo applications.

All-in-One Platform

The 7000 Series has the ability to place implants up to 80 Ncm with a Mont Blanc® 20:1 handpiece.

For oral surgery, predictable and adjustable torque adds to practice efficiency.

For endodontics, the 7000 Series features auto-stop-reverse functionality to prevent file separation and comes preloaded with the complete Dentsply Sirona rotary file library plus the ability to customize your own presets.

For restorative procedures, the 7000 Series streamlines workflow; simply change out the handpiece and you have a restorative powerhouse.

Dynamometer Handpiece Calibration

The 7000 Series dynamometer calibration system gives dentists the ability to quickly and easily diagnose the health of a handpiece and verify that the handpiece is producing the correct torque for the application.

Upgradeable Software

Upgradeable software means a fantastic return on investment while ensuring that your future needs are also met.

Autoclavable Tubing

Reusable, autoclavable irrigation tubing means less waste, no more buying or keeping track of costly disposables.

Multi-Function Foot Control

The Multi-Function Foot Control replaces the need to adjust the unit manually during the procedure providing the operator variable speed control, forward-reverse, irrigation flow rate, torque adjustment, and memory functions.

Available handpiece ratios

1:1, 1:2, 1:3, 1:5, 8:1, 20:1

AEU-7000

- Recommended for oral surgery, dental implants, rotary endodontics & restorative dentistry
- All-in-one platform streamlines workflow
- Compatible with E-type 8:1 endo, 20:1 implant, 1:1, 1:2 surgical, 1:5, 1:3 high speed handpieces with external irrigation connections
- Additional handpiece ratios available via dynamometer calibration 4:1, 16:1, 32:1
- Advanced features and efficient adjustable torque desired by specialists
- Dynamometer provides exceptional handpiece testing and calibration
- Autoclavable 40k rpm brushless micromotor
- Works for all implant types
- Endodontic presets and complete Dentsply Sirona rotary endodontic file library onboard
- Upgradeable software
- Adjustable irrigation flow rate of up to 140 ml/min.
- Reusable, autoclavable irrigation tubing
- Easy cleaning and maintenance
- Optional: fiber-optic illumination, multifunction foot control, carrying case
- Console dimensions 9.98" W x 9.42" L x 5.1" H (25 x 24 x 13 cm)
- Universal voltage 100-240V

7000 Motor with On/Off Foot Control	AEU-7000E	\$5,230
7000 Motor with Multi-Function Foot Control	AEU-7000E-70V	\$5,450
7000 Motor with LED Illumination and Multi-Function Foot Control	AEU-7000L-70V	\$5,750

Options

Carrying Case	410205-HNDL	\$950
Micromotor with 10ft cable	AE-230M-40EXT	\$1,200
Fiber Optic Micromotor with 10ft cable*	AE-230L-40EXT	\$1,420
On/Off Foot Control	AE-7PM	\$190
On/Off Foot Control with 10ft cable	AE-7PM-10	\$210
Multi-Function Foot Control with 10ft cable	AE-70V2-10	\$660

Replacement Parts

Fiber Optic Micromotor*	AE-230L-40	\$1,300
Micromotor with 6ft cable	AE-230M-40	\$1,100
Multi-Function Foot Control	AE-70V2	\$610
Calibration Bur	461558	\$16
Irrigation Bag Hanger Rod	461541	\$48
Power Cord	840079	\$23
Handpiece Cradle	461561-01	\$11
Irrigation Tubing Set	AE-23	\$80

*Requires AEU-7000L Console



"With such a heavy motor workload daily it was imperative to pick something with enough cutting power for both procedures daily. The Aseptico 7000 gave us some unique features so this is one piece of equipment that I would not want to practice without."

Dr. Paresh Patel



Watch a Video
Demonstration



Blog: Aseptico® 7000
Series Dental Motor is
Named Top 25 Implant



PACKAGE DISCOUNTS AVAILABLE.
CALL TODAY! 425.408.9347

6000 Series Motor All-In-One Platform Our Best Value!

Osteotomy & Implant Placement,
Rotary Endo & Restorative

AEU-6000

- Recommended for implantology and rotary endodontics
- All-in-one platform streamlines workflow
- Compatible with E-type 8:1 endo, 20:1 implant, 1:1, 1:2 surgical, 1:5, 1:3 high speed handpieces with external irrigation connections
- Integrated Handpiece Calibration
- Autoclavable 40k rpm brushless micromotor
- Works for all implant types
- Endodontic presets
- Upgradeable software
- Adjustable irrigation flow rate of up to 140 ml/min.
- Reusable, autoclavable irrigation tubing.
- Easy cleaning and maintenance
- Optional Multi-Function Foot Control
- Small footprint:
8.6 W x 9.2" L x 4.8" H (22 x 23 x 12 cm)
- Universal voltage: 100-240V

The AEU-6000 is our best-selling all-in-one motor because of its reliability and affordability. This versatile machine has the ability to place implants up to 50 Ncm with a Mont Blanc 20:1 handpiece. It performs rotary endo with auto-stop-reverse, and it can also be used as a backup for your restorative and oral surgery systems. No other motor on the market can give you these benefits at such an affordable price!

With upgradeable software, the AEU-6000

unit will be a workhorse for many years to come. Autoclavable irrigation tubing means you don't have to buy a new set for each procedure. Five preset buttons can be programmed separately for both implants and endodontics, reducing your procedure time dramatically.



"I have really enjoyed my Aseptico AEU-6000. It meets all my surgery needs from placing implants to removing full bony impacted wisdom teeth. I have found that this unit has more features and better torque than my old W&H."



Heath Snell, DDS
General Dentist

Implant Motor w/ Multi-Function Foot Control	AEU-6000-70V	\$3,950
Implant Motor w/ On-Off Foot Switch	AEU-6000	\$3,730

Options

Carrying Case	410205-HNDL	\$950
Micromotor with 10ft cable	AE-230M-40EXT	\$1,200
On/Off Foot Control	AE-7PM	\$190
On/Off Foot Control with 10ft cable	AE-7PM-10	\$210
Mult-Function Foot Control with 10ft cable	AE-70V2-10	\$660


Replacement Parts

Micromotor with 6ft cable	AE-230M-40	\$1,100
Multi-Function Foot Control	AE-70V2	\$610
Irrigation Bag Hanger Rod	461541	\$48
Power Cord	840079	\$23
Handpiece Cradle	461561-01	\$11

Tubing and Accessories see page 20-21

Electrical Operative Motor Top-Rated Motor Converts Air to Electric

AEU-5000

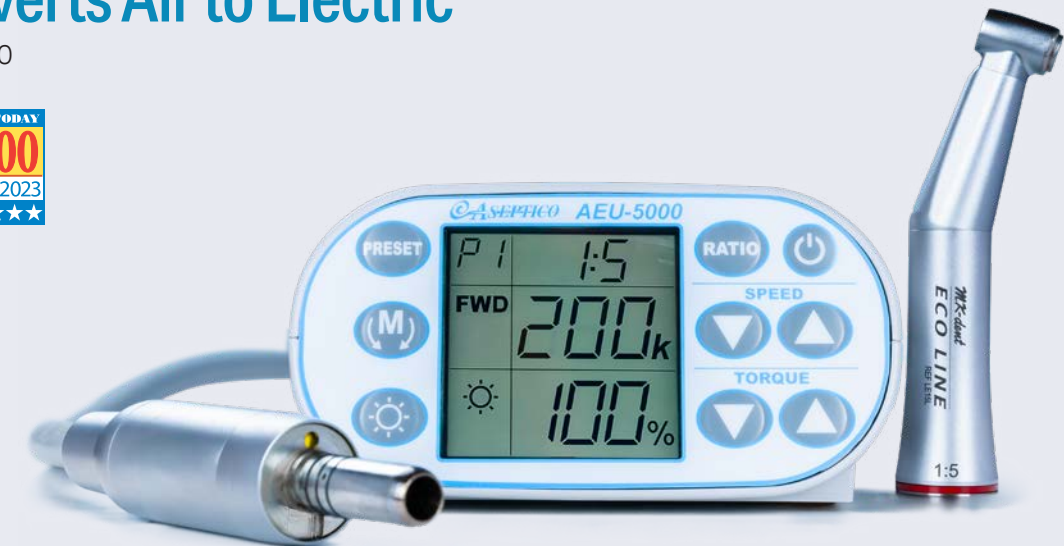

Recommended handpiece 1:5 Restorative with fiber-optic illumination (AHP-72MK-FO)

- Converts air to electric
- Powerful, high torque electric motor
- Attach to your air-driven high speed connection and use your rheostat
- Electrical option for efficient and quiet restorative procedures
- Up to 200,000 rpm with high speed handpiece
- Intuitive design and operation
- Includes endo mode with auto-stop-reverse torque setting
- Universal voltage 100-240V

The versatile and cost-effective operative motor

is designed to convert a standard air dental delivery system into a powerful electric handpiece system for numerous dental applications.

The AEU-5000 features endo functionality, fiber optic intensity adjustment, a compact and easy-to-clean console, a large backlit display, intuitive user interface, and a 40,000 rpm brushless autoclavable micromotor. Innovative LED technology produces the brightest illumination you'll find in a dental handpiece.

Electrical Operative Motor **AEU-5000** **\$1,890**

Adjustable Mounting Bracket **330569** \$70

Replacement Micromotor **AE-240-40** \$1,450

Recommended Handpieces see pages 14- 18



Available handpiece ratios 1:5, 1:2, 1:1, 5:1, 8:1



FRIENDLY AND PROFESSIONAL CUSTOMER SERVICE. CALL TODAY! 425.408.9347

CALL TODAY! 425.408.9347 or shop online at aseptico.com



Endo DTC® Rotary Motor America's Most Popular

AEU-25

- Most sold, most trusted, endo motor in America
- Dentsply Sirona file library included
- Auto-stop-reverse functionality



Foot control included



"As a busy endodontist, I want equipment that is reliable, but more importantly, from a company that has great customer service. Aseptico delivers on both fronts. Fortunately, their handpieces are workhorses, so I don't have to contact them very often!"

Chad Burgess, DDS
Endodontist

The Endo DTC Rotary Endo Motor is America's most popular endo motor because it is a proven reliable performer. Many dentists did their endodontic training on this system and they still prefer it today. And why not? There is so much packed into this small footprint at such an affordable price: advanced digital torque control technology with auto-stop-reverse helps prevent file separation, five programmable memory buttons, onboard file library, integrated calibration system, autoclavable 30k rpm brushless micromotor, and upgradeable software. Pair with a MicroNiTi series handpiece to make the perfect rotary endo system.

Available handpiece ratios
1:5, 1:1, 4:1, 8:1, 16:1, 64:1

Endo DTC®	AEU-25	\$2,200
Replacement Parts		
Micromotor	AE-4B-30	\$770
On/Off Foot Control	AE-7P	\$190
Handpiece Holder	AE-5FG	\$24

Recommended handpieces see page 15



MicroNiTi Handpieces

We recommend the MicroNiTi series of endodontic handpieces. AHP-88MN & AHP88-MNP). Select from Minihead or Microhead. See page 15.

www.aseptico.com



Electronic Endodontic System Our Most Advanced Endo Motor

AEU-27A

- Lightweight motor and cable
- Upgradeable software
- Weighted omnidirectional foot pedal
- Includes all features of America's most popular endo motor, the Endo DTC®



The Electronic Endodontic System features digital torque control with auto-stop-reverse to prevent file separation and automatic calibration to test both the motor and handpiece for accurate performance and consistent operation. Easy to operate, with an auto-clavable brushless motor, this motor also features onboard file library, five programmable memory locations, and preset options that allow single file, file series, or custom file series. Upgradeable software assures return on investment. Includes weighted omnidirectional foot pedal.

Available handpiece ratios
1:5, 1:1, 8:1, 16:1

Electronic Endodontic System	AEU-27A	\$2,300
Replacement Parts		
Micromotor	AE-250-30	\$670
On/Off Foot Switch	AE-7PM	\$190

MicroNiTi series endodontic handpieces recommended see page 15

Richard Mounce, DDS
Mounce Endo
Eagle River, AK

"The AEU-27A is the most versatile, safe and economical motor available anywhere in the world, at any price. That's why I customized the AEU-27A with rpm and torque presets for the MounceFiles in both Controlled Memory and Standard Nickel Titanium. Upgradeable, simple to use, and ergonomic, it sets the standard for electric endodontic motors for many years to come."



Don't forget the HandiDam®

Created by Aseptico, it's the original pre-framed rubber dam trusted worldwide. Easier and more efficient — it's the only rubber dam that allows access to suction or take x-rays while isolation and clamps remain in place. See page 28.



FRIENDLY AND PROFESSIONAL CUSTOMER SERVICE. CALL TODAY! 425.408.9347

Mont Blanc® World Class Implant Handpieces

The Mont Blanc Implant Handpiece is crafted from stainless steel and titanium, making it lightweight and durable. Precision ball bearing construction reduces friction to prevent overheating during surgical procedures. It features smooth surfaces that improve infection control and ensure easy cleaning, with the ability to remove the head for cleaning without tools.

Electric (E-type) | 80 Ncm max torque | 20:1 reducer | 15-2,000 rpm (40k rpm micromotor)

Impulsion Best Value Implant Handpiece

Rated at 50-55 Ncm, the 20:1 Impulsion Dental Implant Handpiece offers great value. This reliable handpiece uses sleeve bearings, making it a good handpiece for drilling an osteotomy at an affordable price.

Electric (E-type) | 55 Ncm max torque | 20:1 reducer | 15-2,000 rpm



AseptiSpray™ Ensuring Handpiece Longevity



Because of our exceptionally high standards, and because proper handpiece maintenance is so important for performance and longevity, Aseptico has been making its own handpiece lubricant for years to make certain that dentists have an excellent formulation for today's dental handpieces. AseptiSpray™ Advanced Handpiece Cleaner and Lubricant is designed for complete internal maintenance of turbine, air motor, contra angle, and E-type handpieces. Aseptico's formulation contains no CFC's and is ozone safe. For proper cleaning and maintenance of handpieces, it is a must. AseptiSpray™ is affordably priced for the dental practice.



ASEPTICO SELECT HANDPIECES ARE EXCEPTIONAL PRODUCTS,
RIGOROUSLY QUALIFIED AND PERFORMANCE TESTED.

See all handpiece accessories, including AseptiSpray™ handpiece cleaner and lubricant, irrigation clips, cannulas, and more on pages 19-20.

1:1 Straight Surgical

AHP-64MKI

- Electric (E-type)
- 1:1 low speed
- 300 to 40,000 rpm (w/ 40k rpm micromotor)



MK-dent 1:1 Straight Electric Surgical Handpiece is ideal for osteoplasty, working dentin, adjusting dentures, sectioning enamel and performing sinus lifts. The spray nozzle ensures safe working conditions due to the optimum cooling of the preparation site. The nozzles are made to provide a clear view of the surgical area as well as provide a pleasant noise level.

Straight Electric Handpiece	AHP-64MKI	\$620
AseptiSpray™ Cleaner and Lubricant	AHP-09	\$39

20° Surgical

AHP-77W

- Electric (E-type)
- 20° bend for optimum comfort
- 1:2 increaser
- 600 to 80,000 rpm (w/ 40k rpm micromotor)
- Requires 70 mm bur length



AHP-77W 20° Electric High Speed Handpiece is the oral surgeon's number one choice for third molar extractions. The angled design prevents fatigue by allowing it to rest on the back of the hand. A scratch-resistant surface and reinforced nose bearings guarantee durability, and the sliding lip seal protects the interior from penetration of dirt particles.

20° Electric Handpiece	AHP-77W	\$2,047
W&H Cleaner and Lubricant	AHP-10 & AHP-10N	\$76.75
Irrigation Clip	AHP-77W-CLP	\$69.75

45° Surgical

AHP-79W

- Electric (E-type)
- 1:2.7 increaser
- 810 to 108,000 rpm (w/ 40k rpm micromotor)



W&H 45° High Speed Surgical Handpiece is ideal for precise tooth separations and apical resections. The 45° head allows for optimum access to the posterior of the mouth. Available in non-fiber optic and self-generating fiber-optic options, this handpiece comes complete with a triple-spray head that guarantees sufficient cooling of the surgical bur and site.

45° Electric Handpiece	AHP-79W-FO	\$2,123
45° Electric Handpiece	AHP-79W	\$1,761
W&H Cleaner and Lubricant	AHP-10 & AHP-10N	\$76.75

1:2 Straight Surgical

AHP-77DS

- Electric (E-type)
- 1:2 increaser
- Up to 80,000 rpm



AHP-77DS Electric High Speed Handpiece is lightweight and compact. The traditional straight design allows for exceptional posterior access. With its 1:2 increaser and rotation speeds up to 80,000 rpm, this handpiece is a great option for bone contouring and sectioning teeth.

1:2	AHP-77DS	\$1,755
AseptiSpray™ Cleaner and Lubricant	AHP-09	\$39

Minihead Endodontic

AHP-88MN

- Electric (E-type)
- Swing latch
- Speed range (w/ 30k rpm micromotor)
- 8:1 — 40 to 3,600 rpm



Anthogyr MicroNiTi Minihead Endodontic Handpiece has a swing latch head and provides unrivaled accuracy, whether it is being used with a rotary or reciprocating motor. The mini-head size allows for excellent access and line of sight during procedures. Made for nickel-titanium endodontic files, the MicroNiTi Minihead is compatible with all ISO rotary endo files.

8:1	AHP-88MN	\$680
AseptiSpray™ Cleaner and Lubricant	AHP-09	\$39

Microhead Endodontic

AHP-88MNP

- Electric (E-type)
- Push button
- Speed range (w/ 30k rpm micromotor)
- 8:1 — 40 to 3,600 rpm



Anthogyr MicroNiTi Microhead Endodontic Handpiece has a push button head and provides the same unrivaled accuracy as other MicroNiTi series handpieces, whether it is being used with a rotary or reciprocating motor. The micro-head size allows for exceptional access and line of sight during procedures. Made for nickel-titanium endodontic files, the MicroNiTi Microhead is compatible with all ISO rotary endo files.

8:1	AHP-88MNP	\$750
AseptiSpray™ Cleaner and Lubricant	AHP-09	\$39

1:3 Restorative

AHP-65TI-FO and AHP-65TI

- Electric (E-type)
- External irrigation
- 1:3 increaser
- Optional fiber optics
- Max 120,000 rpm (w/ 40k rpm micromotor)
- Perfect for high speed dentistry with your implant/surgical motor



With an external irrigation port and 1:3 increaser gearing, the AHP-65TI Restorative Electric High Speed Handpiece makes any surgical motor with irrigation into a restorative powerhouse. Ergonomic and light-weight design allows for optimum comfort during your procedure.

1:3 High Speed w/ Fiber Optics	AHP-65TI-FO	\$1,600
1:3 High Speed	AHP-65TI	\$1,100
AseptiSpray™ Cleaner and Lubricant	AHP-09	\$39

1:5 Restorative

AHP-72MK-FO

- * Electric (E-type)
- Internal irrigation
- 1:5 increaser
- Fiber optics
- Max 200,000 rpm (w/ 40k rpm micromotor)
- Perfect for high speed dentistry with your electric operative motor



Provides powerful, quiet, vibration-free operation with class-leading precision, the 1:5 Restorative Handpiece features Long Lasting Gears that are engineered to excel at every angle, and ceramic bearings that stand up to heavy use. Fiber optics optimize visibility with 100% glare-free illumination at a color temperature of 5500 K.

1:5 High Speed	AHP-72MK-FO	\$990
AseptiSpray™ Cleaner and Lubricant	AHP-09	\$39

4-Hole Restorative

AHP-96MK

- Air-driven
- Internal irrigation
- Non-optic
- Max 320k rpm, 24 watts
- Max volume 62 dB
- Recommended for high speed dentistry



Simple to use technology, developed for what you do every day. The 4-hole High Speed Handpiece features a clean, elegant shape and a virtually scratch proof smooth chrome surface. Manufactured using proprietary laser balancing technology, this handpiece delivers 24 watts of torque and features a premium cartridge for enduring precision and prolonged life.

MK Dent 4-hole High Speed Restorative Handpiece	AHP-96MK	\$480
AseptiSpray™ Cleaner and Lubricant	AHP-09	\$39

5-Hole Restorative

AHP-97MK

- Air-driven
- Internal irrigation
- Fiber optic glare-free, color accurate illumination
- Max 320k rpm, 24 watts
- Max volume 62 dB
- Perfect for high speed dentistry



It feels great in your hand and sports a clean modern shape with a chrome surface. Fiber optics deliver 100% glare-free, color accurate illumination. The 5-hole High Speed Restorative Handpiece is equipped with the AeroPerformance Blade premium cartridge. Ceramic bearings stand up to heavy use and autoclaving.

MK Dent 5-hole High Speed Restorative Handpiece	AHP-97MK	\$510
AseptiSpray™ Cleaner and Lubricant	AHP-09	\$39

6-Pin Restorative

AHP-82MK

- Air-driven
- Internal irrigation
- Fiber optic with LED color accurate illumination
- Max 350k rpm, 25 watts
- Max volume 55 dB



Excellent for high speed dentistry in-office or with a mobile pneumatic delivery system. More power, better sound, spectacular control. The MK Dent 6-pin High Speed Restorative Handpiece is built for the pro user who demands best in class. Reverse flow protection minimizes the entry of contaminants into the head. Fiber optics provide color-accurate illumination. Supremely finished with a scratch-resistant titanium surface. Dentists will enjoy working with this exceptional handpiece.

MK Dent 6-pin High Speed Restorative Handpiece	AHP-82MK	\$1,000
MK Dent 6-pin Coupler	AHP-82MK-QD	\$370
AseptiSpray™ Cleaner and Lubricant	AHP-09	\$39

6-Pin Fiber Optic Coupler

AHP-82MK-QD

- Air-driven
- Works with LED bulbs and glass rod fiber optics
- Color accurate



AHP-82MK Handpiece shown with AHP-82-MK-QD Coupler attached.

This 6-Pin Fiber Optic Coupler comes with a non-retraction valve and water volume regulation. Works with both LED and glass rod fiber optics to provide 100% glare free and color accurate light.

6-Pin Fiber Optic Coupler **AHP-82MK-QD** \$370

CALL TODAY!
425.408.9347
Or shop online at aseptico.com

Low Speed Straight

AHP-101

- Electric (E-type)
- Attaches to E-type air motor
- Also pairs with AHP-50MK
- 1:1 Straight Handpiece
- Post for Doriot-style prophy angle attachment
- Recommended for use with Aseptico portable systems and lab motors



The Low Speed Straight Handpiece is made of lightweight aluminum and is designed to work with Aseptico® air-driven portable dental systems. The AHP-50 Low Speed Air Motor, which allows air turbines to accept E-type handpieces, is required to connect this handpiece to Aseptico portable dental systems.

Low Speed Straight Handpiece **AHP-101** \$380

Low Speed Contra Angle

AHP-63MK

- Electric (E-type)
- Attaches to E-type air motor
- Also pairs with AHP-50MK
- Swing Latch
- 1:1 Contra Angle
- Recommended for low speed work on portable equipment or in an office operatory



Low Speed Contra Angle Handpiece is a D-latch style low speed handpiece that can be attached to any E-type air motor. Made by MK Dent

MK Dent 1:1 Swing Latch Contra Angle **AHP-63MK** \$260

AseptiSpray™ Cleaner and Lubricant **AHP-09** \$39

4-Hole to E-Type Air Motor

AHP-50MK

- Air Driven
- Low Speed Air Motor (4-hole)
- Adapts an air system to accept an E-type handpiece



The Low Speed Air Motor accepts E-type handpieces. Recommended for use with Aseptico® air-driven portable dental systems or other in-office applications where you need to adapt an air system to an E-type handpiece.

Low Speed Air Motor **AHP-50MK** \$360



AseptiSpray™ Handpiece Cleaner and Lubricant

AseptiSpray™ is a must for proper cleaning and maintenance of handpieces. It keeps handpieces working efficiently and ensures long life and durability.

AseptiSpray™ w/ Nozzle (8.8 oz can) **AHP-09** \$39
 Nozzle only **AHP-09N** \$12

W&H f1 Service Oil Handpiece Cleaner and Lubricant

Synthetic oil of maximum purity, which has been specially made for maintenance of W&H handpieces. Required by W&H warranty policy.

Handpiece Spray w/ Nozzle **AHP-10** \$76.75 & **AHP-09**



Mont Blanc®/Impulsion Handpiece Accessories

Depth Stop Attachment	AHP-85MB-STP	\$130
Black Irrigation Clip	AHP-85MB-CLP	\$45
Cannula	AHP-85MB-CAN	\$30
Cannula Clip Set (Includes Small Clip, Large Clip, Cannula and Tubing)	AHP-07	\$24
Cannula and Tubing	AHP-04	\$19
Small Clip	AHP-05	\$8
Large Clip	AHP-06	\$8
Y Connector	AE-2	\$9

Autoclavable Black Irrigation Clip for Mont Blanc® Handpiece **AHP-85MB-CLP**



Tubing & Accessories

Autoclavable surgical tubing saves money and prevents waste. Disposable tubing sets cost about \$12 per set compared to less than \$1 per use for autoclavable tubing that lasts for approximately 100 uses. Autoclavable irrigation clips are money-savers too.

Irrigation Tubing Set for IV Bags (For AEU-6000 and 7000 series)	AE-23	\$80
Irrigation Clips (qty: 10)	AHP-03	\$22
Irrigation Pump Tubing Section (qty: 10)	AE-23-10	\$110
Irrigation Tubing Set for IV Bags (For AEU-1000 and 1070 series)	AE-24	\$80
Irrigation Pump Tubing Section (qty: 10) (For AEU-1000 and 1070 series)	AE-24-PUMP	\$26
Motor Cap (For AE-4B-10/30)	461084-08	\$10
Motor Plug (For AE-4B-10/30)	460183-08	\$11

Foot Controls

Multi-Function Foot Control

AE-70V2

The Multi-Function Foot Control replaces the need to adjust the unit manually during procedures, providing the operator with variable speed control or on/off pedal, forward-reverse, torque adjustment, irrigation flow rate, and memory functions for implants and endodontics. The Multi-Function Foot Control is available with a standard 6 ft. cable (AE-70V2) or a 10 ft. cable (AE-70V2-10). Made in USA.

Multi-Function Foot Control with 6' cable	AE-70V2	\$610
Multi-Function Foot Control with 10' cable	AE-70V2-10	\$660

On/Off Foot Switch



AE-7PM

On/Off Foot Switch for Aseptico dental motors:
AEU-7000 Series, AEU-6000 Series, AEU-27A, and portable systems AEU-350, AEU-525/AEU-525S. Made in USA.

Foot Switch with 6' cable	AE-7PM	\$190
Foot Switch with 10' cable	AE-7PM-10	\$210

Replacement Micromotors

AEU-7000L micromotor and 6' fiber optic cable	AE-230L-40	\$1,300
10 ft extended cable option	AE-230L-40EXT	\$1,420
6000 Series and 7000 Series Implant Motors micromotor and 6' cable	AE-230M-40	\$1,100
10 ft extended cable option	AE-230M-40EXT	\$1,200
AEU-27A Advanced Endodontic Motor micromotor	AE-250-30	\$670
AEU-5000 Electric Operative Motor micromotor	AE-240-40	\$1,450
Micromotor and Cable for:	AE-4B-30	\$770
AEU-25 Endo DTC®		
AEU-707AV2 Implant/Surgical Motor		
AEU-17BV2 SurgiMotor II		
AEU-10SS Light Duty Lab Motor replacement micromotor	AM-103B	\$220
AEU-1000/1070 Micromotor	AE-230A-40	\$1,100
AEU-1000/1070 Micromotor with 10 ft Cable	AE-230A-40EXT	\$1,200

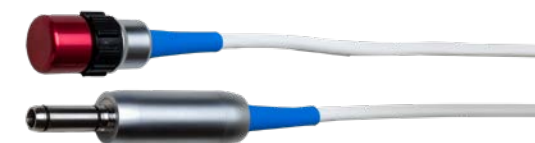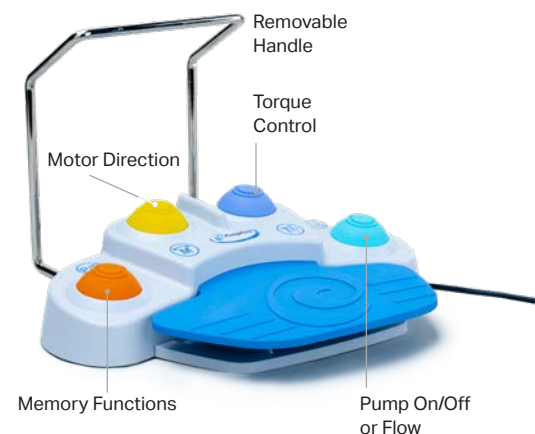

Aseptico

INNOVATIONS
AND UTILITY





Newest Generation Penguin II

Implant Stability Monitor

RFA-55250-MULTI

New and Better! Aseptico is proud to announce the US arrival of the next generation RFA implant stability monitor – the Penguin II.

- Evolved design brings increased comfort, surer grip
- Relocated, easier-to-read LED display screen
- Battery can now be replaced, extending product life
- New vertical wireless charging stand saves space

Implant Stability Monitors

- Provides accurate Implant Stability Quotient measurement for greater confidence determining when to load dental implants
- Cordless, lightweight, and simple to use
- Half the price of competing ISQ devices
- Two-year warranty, Made in Sweden

Penguin Implant Stability Monitors

Now every dentist can access affordable Resonance Frequency Analysis (RFA) technology with the Penguin RFA and Penguin II for assessing osseointegration of dental implants. These devices measure local bone density, translating it into an implant stability quotient (ISQ) value between 1 and 99. Higher ISQ values indicate greater stability. Especially useful for all-on-4 or all-on-X cases, the Penguin monitors provide precise, objective insights into osseointegration, ensuring predictable outcomes

Penguin II	RFA-55250-MULTI	\$2,490
Penguin RFA	RFA-55002-MULTI	\$2,290
Penguin MultiPeg	RFA-550**	\$39

**Part number varies by implant platform.

MulTiPeg™ Sensing Devices

Both Penguin RFA and Penguin II utilize the same MulTiPegs™

- For all major implant systems
- Made of durable, tissue friendly titanium
- Autoclavable for 20 uses
- Laser marked with type number
- ISQ Standard Calibrated™
- Cost effective

Download Entire
MulTiPeg List



Watch a Video
Demonstration



NEW GENERATION

Monitor Osseointegration
Reduce treatment time
Manage risk patients

penguin II



Reduce the Risk of Implant Failure and determine safe loading times with the added confidence of RFA technology. The Penguin II implant stability monitor gives dentists the ability to reliably monitor osseointegration of dental implants.

IMPROVE YOUR IMPLANT PRACTICE
with a **PENGUIN IMPLANT STABILITY MONITOR**
Learn more! Call 425-408-9347 or shop online at aseptico.com

Aseptico

GLOBAL DENTAL CATALOG PROVEN PRODUCTS

Provider Spotlight

MOBILE DENTAL XPRESS

"The GO Gets You Ready to Go and I Love it!" Why Portable Hygienist, Sonya Dunbar, Chooses The GO For Portable Dental Care

Who is Sonya Dunbar what is Mobile Dental Xpress?

Sonya Dunbar, MPH, RDH has been serving in the dental industry for 30 years. With her husband, Gerald Dunbar, they decided to found Mobile Dental Xpress. Their company specializes in bringing dental care to underserved seniors in long term care facilities, giving them the dental attention they desperately need.

"I've been in dentistry for 30 years," says Sonya, "but by far, working with our seniors has been the most rewarding, helping someone that cannot help themselves."

How the GO made an impact on Mobile Dental Xpress

The state-of-the-art GO Ultra Portable Dental System has been a game changer for Sonya and Mobile Dental Xpress. Perhaps the biggest impact has been the small size of the GO and the reduced physical stress it has put on Sonya. "Portable equipment for years has been 60-70 lbs, moving that [older] portable unit out and all your other equipment wears and tears on your body," says Sonya.

At only 32 lbs, the GO Ultra Portable Dental System has changed the mobile dentistry game and is less than half the weight of Sonya's previous mobile dental unit. Such

a difference in weight has allowed Sonya to extend her career due to the reduced stress on her body. This smaller size does not come at the cost of functionality though. Sonya and Mobile Dental Xpress lead dentist, Dr. Marilyn Sandor, DDS, MS, need to regularly be able to provide a full range of dental options for their patients. "Everything you need for that day is there already with the unit, my hygiene equipment, the restorative equip-

ment, even the extraction equipment is there for the dentist," says Sonya. "We love using the GO, we use it whenever we are providing any operative dental services, anything requiring suction and anything that's going to require a handpiece," says Dr. Sandor. The GO Ultra Portable Dental System gives Mobile Dental Xpress increased portability without having to sacrifice care options.



The GO is built to overcome challenges in unique fields

The other huge benefit the GO Ultra Portable Dental System provides Sonya and her team is a reduced noise level. Working in the senior care field, Sonya regularly cares for patients with hearing aids that can be sensitive to loud, motorized equipment. "The GO comes in, it's not

noisy, we're able to be soft, compassionate and provide that dental care they need in their own homes," says Sonya. Crucially, the GO is regulatory compliant. The GO Ultra Portable is made in the USA, FDA-registered, approved by Health Canada, UL-certified, and safe to use in any dental environment.



"I appreciate Aseptico for hearing us, for hearing the mobile clinician, for hearing the public health hygienists and dentists, and making something that can go with our lifestyle,"

**Have Questions? Call us at 425-408-9347
or visit [Aseptico.com/GO](https://www.aseptico.com/GO)**



What is the Industry Saying?

ASEPTICO PRODUCT RECOGNITION



DENTAL PRODUCTS REPORT
TOP 5 IMPLANT PRODUCTS

Aseptico® 1070 Series Implant Motor

Dental Products Report magazine named the Aseptico 1070 Series Implant Motor as one of the **Top 5 Implant Products of 2023!**



For Value & Excellence GoAseptico!

Aseptico® Dental Motors Are Made in America and Trusted Worldwide.

For a Complete Listing of Aseptico Industry Awards, visit aseptico.com



Aseptico® Products Have Been Recognized Multiple Times

by leading dental publications for their innovation and preference for use by dental professionals. *Dentistry Today* and *Dental Products Report*, leading US magazines for dentists, recently recognized multiple Aseptico® products.



Aseptico® 7000 Series All-in-One Motor

- Reader's Choice Top 25 Implant Product March 2023
- Reader's Choice Top 25 Endodontic Product April 2023

Top 25 Implant Product and Top 25 Endodontic Product for multiple years running!



AMC-15 Mobile Dental Cart

- Reader's Choice Top 50 Technology Product Nov/Dec 2023
- Reader's Choice Top Innovative Products Jan/Feb 2023

Spotter Smart Implant Detector Find Buried Implants Faster and Preserve Soft Tissue

SPT-5111001-KIT

- Painlessly and accurately locate implants covered by gum tissue!
- Non-invasive, no exploratory flap needed
- Locate buried implants quickly and easily
- Minimize the recovery time for the patient
- Two-year warranty
- Works with any titanium implant

Spotter Smart Implant Detector locates dental implants covered by gum tissue quickly and painlessly — no exploratory flap needed! Just move the wand over the site and the Spotter will guide you to the exact center of the implant.

Instead of opening a flap to locate the implant, the Spotter finds the exact center quickly, allowing you to punch and reveal only the implant thereby helping to preserve soft tissue and minimize recovery time for the patient.

The Spotter is especially helpful in the case of multiple implants when implant sites may be hard to visualize. The Spotter Smart Implant Detector works with any type of titanium implant.

It is easy to use, requiring only a few minutes to learn how to operate. The unit operates on one AA battery and its detachable sensing head is autoclavable.

Spotter SPT-5111001-KIT \$700
Sensors 5/PKG SPT-5111101 \$130



"With Spotter, we don't have to raise a flap, we can go in a flapless procedure, just punching the soft tissue to get to the center of the implant."



Ziv Mazor, DMD



CALL TODAY! 425.408.9347 or shop online at aseptico.com

HandiDam® pre-framed dental dam




- Pre-framed, flexible design facilitates access to the oral cavity for suction or X-rays
- Superior elasticity and tear-resistance
- Saves time—eliminates the need to remove and replace traditional dam during the procedure
- Ergonomic and comfortable for greater patient acceptance
- Meets FDA Powder-Free standard
- Extremely low protein content reduces patient irritation (<50 micrograms)
- Available in Latex or Latex-Free

HandiDam® Latex (Blue)

Qty 20,	I-70BM	\$41
Qty 100,	I-70BM-100	\$160

HandiDam® Latex-Free (Mauve)

Qty 20,	I-70LF	\$65
Qty 100,	I-70LF-100	\$270

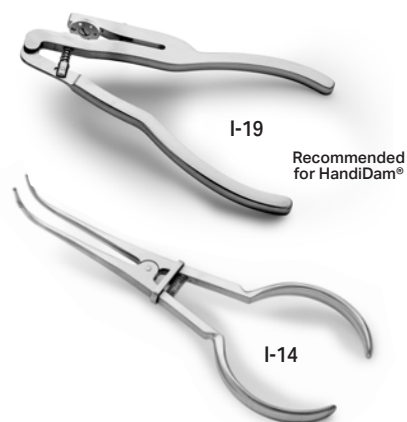

Dental Dam Punch

Ivory Style Punch	I-19	\$230
-------------------	-------------	-------

Dental Dam Forceps

University of Washington Style Forceps	I-14	\$150
--	-------------	-------


TECHNIQUE CHAIN & GRID



Technique Chain & Grid easily attaches to a pre-emerged tooth, so that an orthodontic device may be used to draw the tooth into proper position. Aseptico is the originator of this remarkable traction system that professionals use worldwide.

- 10 sets per box
- 14k gold chain
- 3 grid shapes available

Cuspid (qty: 10)	AMA-09A	\$350
Bicuspid (qty: 10)	AMA-09A2	\$350
Incisor (qty: 10)	AMA-09A3	\$350

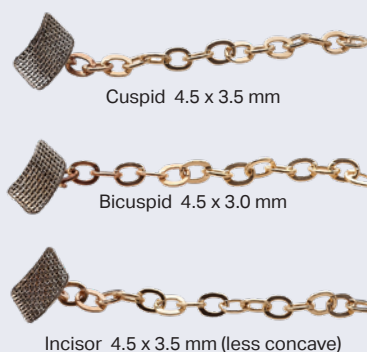

Dental Dam Clamps

Our five most popular clamps.
See full selection at aseptico.com
Dental Dam Clamp \$16 each



K-211
Universal labial cavities on anteriors



K-5
For lower molars



K-9
Most popular anterior clamps



K-90N
Most popular anterior clamps



K-8A
For lower molars, with strong festooned beaks



Compricart II Delivery System

ADU-20B

- One high speed and one low speed handpiece connection
- High volume evacuator and saliva ejector
- Large stainless steel instrument tray (13.5" x 9.75")
- All day water system with two liter self-contained water bottle

Roll the Compricart II Dental Cart from room to room with ease! Great for hospitals and expanding dental offices, the simple controls of the Compricart II make operation stress-free. A compact design, ample work surface, stainless steel instrument tray, integrated water supply, and venturi suction (high volume and saliva ejector) round out the features of this versatile platform. Just connect to an air source – no central vacuum required! Including the optional AWS-1C allows an external scaler to use the cart's water supply.



Compricart II	ADU-20B	\$3,600
Water Connection for External Scalers	AWS-1C	\$85

Options

Fiber Optic Upgrade Kit 5-hole	AA-19LED-04	\$690
Fiber Optic Upgrade Kit 6-pin	AA-19A-06	\$690

Replacement Parts

Complete HVE Assembly	330461	\$520
Disc Foot Control w/ Water Toggle	AA-43W	\$110
Air/Water Syringe	TA-90D	\$85
2L Water Bottle	730665	\$14



AEU-5000 Add-on Option Available
See details on page 9

FRIENDLY AND PROFESSIONAL SERVICE. CALL TODAY! 425.408.9347





Automated Air Delivery System

ADU-07 and ADU-07A

- Choose 2 or 3 handpiece system
- Includes 3-way air/water syringe
- Wall or cabinet mount
- Optional fiber-optics
- Includes foot control


The Automated Air Dental Delivery System is a pneumatic handpiece control system that gives you everything you need at a price you can afford. Highly customizable, select from a 2-handpiece system with manual control, or 3-handpiece system with automatic selection and individual air and water adjustments. Both models include air pressure gauge, 3-way air/water syringe, foot control, and water flow adjustments. Flat top design is ideal for mounting into cabinets. Requires connection to air and water sources.

Options

2-Handpiece/Manual Control	ADU-07	\$1,050
3-Handpiece/Automatic	ADU-07A	\$1,150
5-hole Fiber Optic Upgrade	AA-19LED-04	\$690
6-pin Fiber Optic Upgrade	AA-19A-06	\$690

Replacement Parts

Disc Foot Control w/ Water Toggle	AA-43W	\$110
Air/Water Syringe	TA-90D	\$85
Autoclavable Tips for Air/Water Syringe	TA-1	\$9

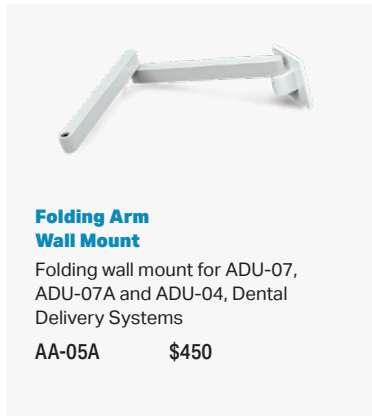

Folding Arm Wall Mount

Folding wall mount for ADU-07, ADU-07A and ADU-04, Dental Delivery Systems

AA-05A \$450



Manual Air Delivery System

ADU-04

- 2-handpiece control system
- 3-way air/water syringe
- Optional independent water
- Ideal for arm or wall mounting
- Requires connection to air

The Manual Air Delivery System is a pneumatic handpiece control system for two handpieces. It includes a 3-way air/water syringe and disk foot control. Requires connection to air and water sources or use optional independent water.

Manual Air Control System	ADU-04	\$620
Clean Water System with 1 Liter Bottle (pictured)	NWS-71	\$130

Replacement Parts

Disc Foot Control	AA-42	\$65
Air/Water Syringe	TA-90D	\$85
1L Water Bottle (for NWS-71)	730471	\$13
Autoclavable Tips for Air/Water Syringe	TA-1	\$9

Powered Trolleys

ATC-12V2 and ATC-15

- Keeps equipment organized and solves operator space concerns
- Storage bin for foot pedals or other items
- Dual outlet with coiled AC power cord
- Adjustable height to easily store under a counter: 26.4" to 37.4" (67-95 cm)
- Four shelves w/ 3 removable stainless steel trays

Aseptico Powered Trolleys give dental teams the ability to easily roll dental motors, scalars, lasers, surgical systems or other electrical equipment from operator to operator. Available in large and standard sizes, each Powered Trolley is adjustable in height and features four shelves of storage with removable stainless-steel trays that fit securely into each shelf. High-quality, sturdy trolleys move and turn with ease and fit into almost any small space.

Standard Powered Trolley ATC-12V2 \$795

Cart size: 26.4-37.4" x 15.75" x 16.5" (67-95 x 40 x 42 cm)
Tray size: 11.3" x 9" (28.8 x 22.8 cm)

Large Powered Trolley ATC-15 \$950

Cart size: 26.5-38" x 19.5" x 19.75" (68-97 x 50 x 50 cm)
Large Tray: 15.4" x 11" (39.1 x 28 cm)



Aseptico Powered Trolleys feature two heavy duty electrical outlets that are certified to safety standards and a sturdy storage bin designed to hold foot pedal or other items.



Our Newest Cart Deluxe Operator Cart

ATC-16

Efficiency Enhancements

The cart's organizational design and mobility reduce clutter and inefficiencies, allowing dental professionals to focus more on patient care.

Hygiene and Maintenance Ease

With its smooth metal and glass surfaces, the cart supports stringent hygiene protocols, simplifying cleaning and sterilization processes.

Aesthetic Appeal

The cart's modern, sleek design elevates the professional appearance of your operator, positively influencing patient perceptions.



Deluxe Operator Cart **ATC-16** \$4,295

Options

Adjustable, Articulating Arm	AA-16ARM	\$600
Motor and IV Bag Holder	AA-16-IMPACC	\$500
Piezo Handpiece Holder	AA-16-PZACC	\$370



270° Safety Glass Door opening with soft-close technology

2 pull out trays



Galvanized steel with stainless steel accents

Ball-bearing casters for easy and smooth movement

Ample tabletop space for all of your equipment



Integrated outlets and automatic cable winding for cable management



3 Drawers with soft close technology



MOBILE & PORTABLE DENTAL EQUIPMENT



2 YEAR WARRANTY
on Aseptico® Portable Equipment Made in the USA*

*1 Year warranty on handpieces. Read our warranty policy on page 33



DENTISTRY ANYWHERE



GO Ultra Portable Dental Delivery System

Our Most Advanced Technology

AEU-350

- Simple to operate for restorative & hygiene
- Integrated brushless motor with LED illumination
- Powerful vacuum and air flow
- Exceptionally quiet operation
- Ultra-compact and lightweight
- Meets carry-on luggage requirements
- Air filter reduces airborne particulates
- Use anywhere in the world with electric or battery power
- Customer-replaceable bacterial-viral filter
- **NEW** Health Canada Approved

Designed for Ease-of-Use and Reliability

The GO leverages Aseptico's expertise in crafting leading portable dental units for today's challenges and opportunities. Designed for operator ease-of-use and outstanding reliability, the GO is capable of restorative dentistry and hygiene functions.

At only 32 pounds, the GO system is fully ten pounds lighter than the Aseptico® Express system — and so compact it qualifies as carry-on luggage for air travel. The GO delivers powerful airflow and suction yet operates more quietly and cleanly than standard air compressor units (it even has an air filter to reduce airborne particulates).

It can operate anywhere in the world — and if there is no electrical power, there is an option for it to operate on a 24-volt battery.

GO Ultra Portable Dental Delivery System **AEU-350** **\$7,650**

Options

24 VDC Battery **AA-350-BATCAB2** \$110
 Power Cord

Replacement Parts

Micromotor w/5'cable	AE-242SC-40	\$1,260
Water Supply Bottle	AA-350BOT	\$22
Waste Collection Tank	330811	\$95
Air Filter	AA-350-FIL	\$37
HVE Valve and Hose Assembly	AA-350-HVE	\$160
Saliva Ejector Valve and Hose Assembly	AA-350-SAL	\$60
Foot Switch	AE-7PM	\$190
Power Cord	840041	\$20
Power Cord Pouch	410222	\$14
Hose Storage Pouch	410225	\$2

Powerful Electric Motor

5k to 200k rpm range is ideal for prophylaxis and restorative.

Faster Drying for Sealants

- Syringe Performance:
- 80 psi initial
 - Up to 50 psi continuous

More Powerful Suction!

New vacuum blower technology

Quieter Operation

A quieter vacuum system that turns on only when HVE and saliva ejector are removed from their holders creates a quieter experience for both operator and patient.

1 Liter Waste Capacity

Easy, clean removal of waste

500 mL Water Capacity

Plenty of water capacity for restorative cases

Impact Resistant

Durable, wheeled case with comfortable handles

Weighted & Omnidirectional Foot Pedal

Stays in place and can be activated from any angle



Super Compact!

22" x 13" x 9.2" fits into the airline overhead as carry on luggage



Easy and Light

Only **32 pounds!** 13 pounds lighter than the Aseptico® Express, the world's #1 system!

Plug in Anywhere

100-240 VAC & 24VDC input options

Can Operate on External Battery Power

Powering the GO Ultra Portable Dental System from a 24 VDC battery is an option for those who work off-the-grid. The cable is purchased from Aseptico, and the customer supplies the 24V battery. The GO can operate for as long as the battery type and capacity will allow.



Sonya Dunbar and The GO

Geriatric Toothfairy | Mobile Dental Xpress



Watch Sonya Dunbar's Testimonial on the GO



OUR MOST ADVANCED TECHNOLOGY!
 CALL TODAY! 425.408.9347

Express Portable Dental System

ADU-17A and ADU-17A2

ADU-17X and ADU-17X2

- Better and Faster Sealants
- Available with one or two air handpiece connections
- High volume evacuator and saliva ejector
- Built-in maintenance free compressor
- Water port for Cavitron/scaler
- High capacity 1 gallon waste container
- 110V or 230V switchable voltage
- 44-48 lbs (20-21.8 kg)
- 10.5" x 14.5" x 24" (27 x 36 x 61 cm)

Aseptico's most popular portable dental systems, the Express Air and Express II portable dental systems are fully self-contained, lightweight, reliable performers that provide everything a dentist or hygienist needs to perform effectively in the field. These powerful units feature up to two air-driven handpiece connections, high volume and saliva ejector suction, internal water reservoir, 3-way syringe, and integrated air compressor. They are transported using a pullout handle and built-in wheel assembly. The system's compact size allows for convenient storage, while the impact-resistant case makes the Express ideal for any field environment. Designed for global use, all Express systems are compatible with 110V or 220V power sources and come in under the airline baggage weight limit.

Express Air – features variable rheostat foot pedal. The handpiece only operates when the variable rheostat foot pedal is depressed; handpiece operates independently from suction and 3-way syringe. **(Dentist preferred)**

Express Air one handpiece connection	ADU-17A	\$5,650
Express Air two handpiece connections	ADU-17A2	\$5,750

Express II – features on/off electric foot switch. All components are activated by the on/off foot switch. This configuration is sometimes desirable for special needs patients or other specialized applications. Please call for more information. **(Hygienist preferred)**

Express II one handpiece connection	ADU-17X	\$5,650
Express II two handpiece connections	ADU-17X2	\$5,750

Fiber Optic Upgrade Kits


5-hole	AA-19LED-04	\$690
6-pin	AA-19A-06	\$690
6-pin LED handpieces	AA-19EGEN-04T6*	\$260

*Will not function with incandescent bulbs.
Consult your Aseptico sales rep for details.



Handpieces sold separately
Recommended handpieces see pages 16-18.

"We love the ease of this piece of Aseptico equipment! It's essential for our business to run efficiently being able to utilize everything we would have in a private office setting while providing care in long-term care facilities. No sacrifice in power or suction either!"



Liz LaJoie, RDH
Director of Clinical Field Teams
Aria Care Partners

Transport III Portable Dental System

AEU-525S and AEU-525

- High torque electric motor with restorative and endodontic programs
- Fiber optic electric handpiece connection
- High volume evacuator and saliva ejector
- Quiet, intuitive, powerful
- Built-in maintenance free compressor
- Water port for Cavitron/scaler
- 115V or 230V Dual voltage
- Optional fiber optic piezo scaler
- 57-63 lbs (25.8-28.6 kg)
- 17.9" x 10.4" x 22.1" (45.5 x 26.4 x 56.1 cm)

The Transport III is Aseptico's most capable portable dental system. Features include noise reduction technology, an intuitive user interface, and 40k rpm brushless auto-clavable micromotor. It is easy to set up, operate, and tear down. All components fit neatly into the impact resistant carrying case with inset handle and wheels.

Delivery module and user interface panel mount externally, placing controls and instruments at a comfortable level.

Dual level fluid sensors alert the user when the waste tank is 85% full for better procedure planning. The Transport III will then shut off the compressor before the tank exceeds capacity. The externally mounted waste container and optional amalgam separator means faster and cleaner waste disposal.

Transport III w/Scaler	AEU-525S	\$12,900
Transport III	AEU-525	\$11,300

Options

Amalgam Separator Kit	AA-525AMAL	\$720
Micromotor w/ 6' cable	AE-240SC-40	\$1,450
Foot Switch with 10ft cable	AE-7PM-10	\$210

Replacement Parts


Micromotor w/5' cable	AE-240SC-40	\$1,450
Waste Tank Assembly	330603	\$660
Foot Switch	AE-7PM	\$190
Power Cord 110V	840049	\$40
Power Cord 220V	840007	\$12
Water Reservoir (white)	730471	\$13

Replacement Scaler Tips For AEU-525S Only

Slim Universal	ASC-10-PE38	\$170
Thin Subgingival	ASC-10-PE37	\$180
Power Universal	ASC-10-PE39	\$130



Handpieces sold separately
Recommended handpieces see pages 16-18.

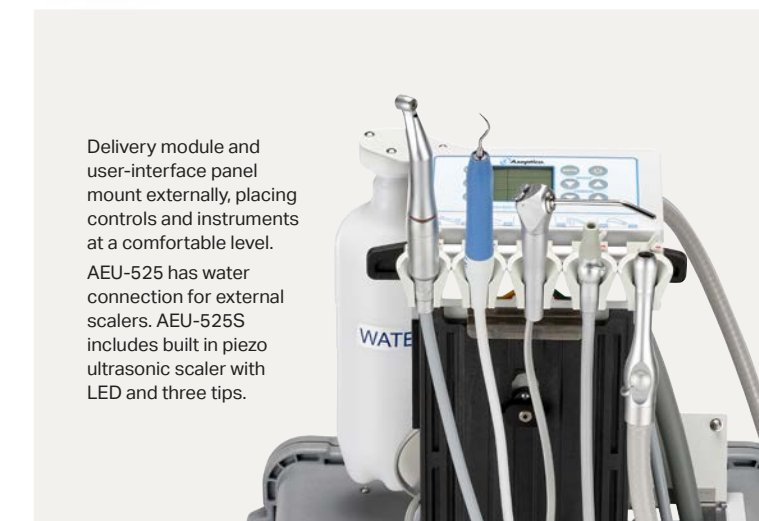

Delivery module and user-interface panel mount externally, placing controls and instruments at a comfortable level. AEU-525 has water connection for external scalers. AEU-525S includes built in piezo ultrasonic scaler with LED and three tips.



ASEPTICO PORTABLE SYSTEMS ARE MADE IN THE USA. CALL TODAY! 425.408.9347

Taskforce Portable Dental System

ADU-12DCE

- New look and case design
- Lightweight and easy to carry, under 33lbs!
- Two air handpiece connections
- Powerful saliva ejector vacuum
- Built-in maintenance free compressor
- Water port for Cavitron/scaler
- 110V or 220V switchable voltage
- 12.5" x 17.5" x 7.5" (31.8 x 44.5 x 19.1 cm)

The most compact pneumatic portable dental system, the Taskforce, represents the best value on the market. At 32.6 pounds the Taskforce is easily carried from location to location. This unit includes two air-driven handpiece connections (high speed /low speed), low volume suction, and air water syringe, making it the ideal system for mobile hygiene. Simple and easy to use, the Taskforce can be set up or put away in just minutes.



Handpieces sold separately

Taskforce **ADU-12DCE** **\$4,150**

Optional Accessories

Handpiece Extension Arm	AA-12A	\$80
Stainless Steel Cupidor with Solids Trap 5" (12.7mm) diameter for convenient expression during treatment	AA-26SS	\$65
High Volume Tip Adapter enables use of standard 11mm high volume tips	AA-12HVE	\$40

Fiber Optic Upgrade Kits

5-hole	AA-19LED-04	\$690
6-pin	AA-19A-06	\$690
Fiber Optic for 6-pin LED handpieces	AA-19EGEN-04T6*	\$260

*Will not function with incandescent bulbs. Consult your Aseptico sales rep for details.

Replacement Parts

Lid Standoff Bracket	AA-12CSB2	\$44
----------------------	------------------	------

"I have used the Aseptico Taskforce Deluxe in 25 developing countries on dental missions. The unit is lightweight, and has always functioned well with minimal power requirements. I have recommended its use to many needing a lightweight, reliable dental operating unit."



COL (Ret) Robert Meyer, DMD
Executive Director
Christian Dental Society



4-Hole to E-Type Air Motor AHP-50
Low Speed Air Motor (4-hole) adapts air turbine to accept E-type handpieces.

Call for more information.



4-Hole Restorative AHP-96MK
Air-driven
Internal irrigation
Non-optic
Max 320k rpm
Max volume 62 dB
Recommended for high speed dentistry
See page 16.

Lightweight Portable Delivery Systems

ADU-04ST and ADU-03ST

- Compact design and case under 15lbs (6.8kg)
- Customizable handpiece and vacuum options
- 3-way air/water syringe
- 25 ft. air line to connect to any air source
- Instrument tray
- Padded carrying case
- Carrying case: 25" x 11.5" x 11.5" (63.5 x 29 x 29cm)
- Requires air compressor for operation
- Quick disconnect options, .25 QD-Standard or 3/8 Schrader (optional)

Restorative System w/ 2 HPC Connections **ADU-04ST** **\$1,890**

Restorative System w/o HPC Connections **ADU-03ST** **\$1,630**

Optional Upgrades

High Volume Suction	AA-21ST	\$730
Saliva Ejector	AA-25ST	\$670
5-hole Fiber Optics	AA-19LED-04	\$690
6-pin Fiber Optics	AA-19A-06	\$690
Fiber Optics for 6-pin LED handpieces	AA-19EGEN-04T6*	\$260

*Will not function with incandescent bulbs. Consult your Aseptico sales rep for details.

The ADU-04ST is designed for dentistry and can run either a high-speed or low speed handpiece (shown with high volume suction and saliva ejector suction options added).



The ADU-03ST is designed for hygiene and sealants (shown with optional high volume suction added).

Airpac III Portable Compressor

AA-100-110

- 1.0 horsepower, 100 psi oil-less motor
- Air delivery 3.1 CFM at 40 PSI
- Built-in wheels and pull handle
- Available in 110V
- Not available in Canada
- Weight: 46 lbs. (20.87 kg)
- Size: 17.93" W x 22.06" L x 10.43" H (45.5 cm W x 56.0 cm L x 26.5 cm H)

The Airpac III Portable Compressor is built by Aseptico to stand-up to heavy use and travel, the compressor is secured in its own super duty carrying case with built-in wheels and pull handle. The Airpac III includes an air storage tank and automatic pressure switch control.

110V Compressor **AA-100-110** **\$2,250**



ASEPTICO PORTABLE SYSTEMS ARE MADE IN THE USA. CALL TODAY! 425.408.9347



AMC-25 Mobile Dental Cart

The **AMC-25 Mobile Cart** is a superior solution for providing general dental care in nursing homes, hospitals, or anywhere a mobile dental unit is needed. Available optional accessories allow for customizing this dental cart to meet your needs, including the ability to add electric handpieces, handpiece illumination, and much more!



Watch CareSouth Carolina's Testimonial on the AMC-25



AMC-25 Mobile Dental Cart

Options	AMC-25	\$14,990
Integrated Electric Motor	330674	\$1,890
Valo Curing Light	ALU-25CLK-VALO	\$1,795
2 Fiber Optic HP Connections		
5-hole (ISO-B)	AA-19LED-25CK5	\$1,320
6-pin (ISO-C)	AA-19A-25CK	\$1,270
Piezo Scaler	ASC-10-25CKL	\$1,750
Amalgam Separator	AA-20AMAL-CK	\$850
Auxiliary Arm with Accessory Tray	AA-20A	\$500
Hospital Kit (external air, water and vac connections)	AA-25HK	\$1,670

Replacement Parts		
Micromotor	AE-240SC-40	\$1,450
Wet/Dry Foot Control	AA-43W	\$110
Air/Water Syringe	TA-90D	\$85
1L Water Reservoir - Clear	730631-01	\$70
1L Water Bottle - Opaque (each)	730471	\$13
Strainer	730745	\$34
Amalgam Canister	730595	\$160

The AMC-25 comes standard with pneumatic handpiece connections, but can be configured for electric handpieces, or a combination of air-driven and electric.

- No plumbing required!
- 2 automatic high speed air handpiece connections (4-hole; non-fiber optic)
- 1 automatic low speed handpiece connection (4-hole; non-fiber optic)
- High and low volume suction with automatic overflow protection and solids trap
- 3-way air/water syringe
- Quick connect for ultrasonic scaler
- 2 hospital-grade outlets for accessories
- Wet/dry foot control
- Dual voltage 110V or 220V
- Operating pressure 80 psi
- High volume vacuum 7.6 SCFM @ 1.3 in Hg, 5.6 SCFM @ 4 in Hg
- Low volume vacuum 2.2 SCFM @ 1.5 in Hg
- Noise Level 65 dBA or less @ 3'4" (1 meter)
- Weight 167 lb (75.8 kg)
- Dimensions 23.5" W x 30" L x 36.5" H (56.7 x 76.2 x 92.7 cm)
- Optional: integrated electric motor, fiber-optic illumination for handpieces, ultrasonic scaler, curing light, amalgam separator, external vacuum/air/water hospital connections, auxiliary arm with accessory tray



AMC-15 Mobile Dental Cart

Designed for Easy and Affordable Mobility, the AMC-15 Mobile Dental Cart is equipped for general dentistry with two automatic high speed air connections, one low speed handpiece connection, internal vacuum and air, water supply, large waste tank with automatic overflow protection, and amalgam separator.

- No plumbing required!
- 2 automatic high speed air handpiece connections (4-hole; non-fiber optic)
- 1 automatic low speed handpiece connection (4-hole; non-fiber optic)
- High & low volume suction with automatic overflow protection and solids trap
- 3-way air/water syringe
- Quick connect for ultrasonic scaler
- Amalgam Separator
- Wet/dry foot control
- Available in 110V or 220V
- Operating pressure 80 psi
- High volume vacuum 6.0 SCFM @ 0 in Hg, 3.6 SCFM @ 4 in Hg
- Low volume vacuum 1.3 SCFM @ 1.5 in Hg
- Noise Level 63 dBA @ 3'4" (1 meter)
- Weight 139 lb (63 kg)
- Dimensions 23.5"W x 30" L x 36.5" H (56.7 x 76.2 x 92.7 cm)
- Optional: auxiliary arm with accessory tray

AMC-15 Mobile Dental Cart

Options	AMC-15	\$11,990
Auxiliary Arm with Accessory Tray	AA-20A	\$500
Extended Fiber Optic Tubing 5-hole (6 ft)	AA-19LED-04	\$690

Recommended Handpieces

Air-Driven see pages 16-17

Replacement Parts		
3-Way Air/Water Syringe with Tip	TA-90D	\$85
Amalgam Separator Filter	730595	\$160
Clear 1-liter Water Bottles w/Level Indicators	730631-01	\$70
Two 1-liter Water Bottles	730471	\$13
Wet/Dry Foot Control	AA-43W	\$110
Power Cord 115V	840049	\$40
Power Cord 230V	840007	\$12



Equipped for general dentistry, the AMC-15 has been recognized for its innovative technology, however, we are most proud of its affordable price — which includes an amalgam separator!



ASEPTICO PORTABLE SYSTEMS ARE MADE IN THE USA. CALL TODAY! 425.408.9347

NEW FABRIC MATERIALS
Better, Stronger, and More Durable



Chairs and Stools Available in Red, Gray and Blue



All seating includes a padded carrying case



AseptiStool

ADC-08

- Set-up and store in seconds
- Includes padded case – no loose parts to misplace!
- 300 lb weight capacity
- Adjustable seat height: 19" - 22" (48.3-55.9 cm)
- 15 lbs (6.75 kg)
- Padded carrying case: 18" x 22" x 8" (48 x 53 x 23 cm)

Durable, adjustable and foldable, the AseptiStool Portable Dental Stool is the perfect seating solution for any clinician. Its one-piece design sets up or stores in seconds, and it swivels on four casters for easy mobility during treatment. It is height adjustable, lightweight, and folds up for simple storage in the included padded travel case.

AseptiStool	ADC-08	\$970
Replacement Padded Carrying Case	410093	\$180



Hydraulic AseptiStool

ADC-10

- Easy set-up and storage
- Includes padded case
- 300 lb weight capacity
- Hydraulic adjustable height: 18" - 25.5" (45.7 - 64.8 cm)
- 19 lbs (8.6 kg)
- Padded carrying case: 26" x 26" x 11" (56 x 56 x 28 cm)

With hydraulic adjustable height, the Hydraulic AseptiStool gives the mobile dentist the comfort and flexibility of doing dentistry in the office. It is easy to set up, tear down and transport in its included padded case. Even with its sturdy construction and high weight capacity, the stool is lightweight and carrying it into any mobile dental scenario is simple.

Hydraulic AseptiStool	ADC-10	\$1,300
Replacement Padded Carrying Case	410213	\$230



AseptiChair

ADC-01


- Sturdy, high-strength steel construction
- Lightweight, only 30 lbs (13.5 kg) in case!
- 300 lb weight capacity
- Backrest adjusts from fully upright to supine
- Price includes padded case, arm slings, and headrest cushion
- Adjustable seat height: 11" - 22" (11 - 55.9 cm)
- Padded carrying case: 23" x 34" x 8" (58 x 86 x 20.3 cm)

The AseptiChair's sturdy rolled steel construction is lightweight yet strong, making it perfect for portable dentistry. It is adjustable to five positions, with seat height ranging from 11 to 22 inches. The toe board folds down, allowing this portable dental chair to be used for examination. The included carrying case makes mobility and transportation simple. No assembly required.


AseptiChair	ADC-01	\$2,200
--------------------	---------------	----------------

Replacement Parts		
Replacement Padded Carrying Case	410193	\$280
Replacement Leg Adjustment Pin	460445-08	\$17
Endcap Heavy Duty 1"	850026	\$3
Armrest Support Pin	460444-08	\$9

"Having quality products like Aseptico patient/operator chairs gives our providers, patients, and facilities confidence in the care that a mobile alternative practice can provide. The support staff at Aseptico is always willing to help when needed. 100% satisfaction guaranteed."

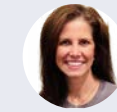

Sara Luetke, RDH
Founder
Sound Dental Care

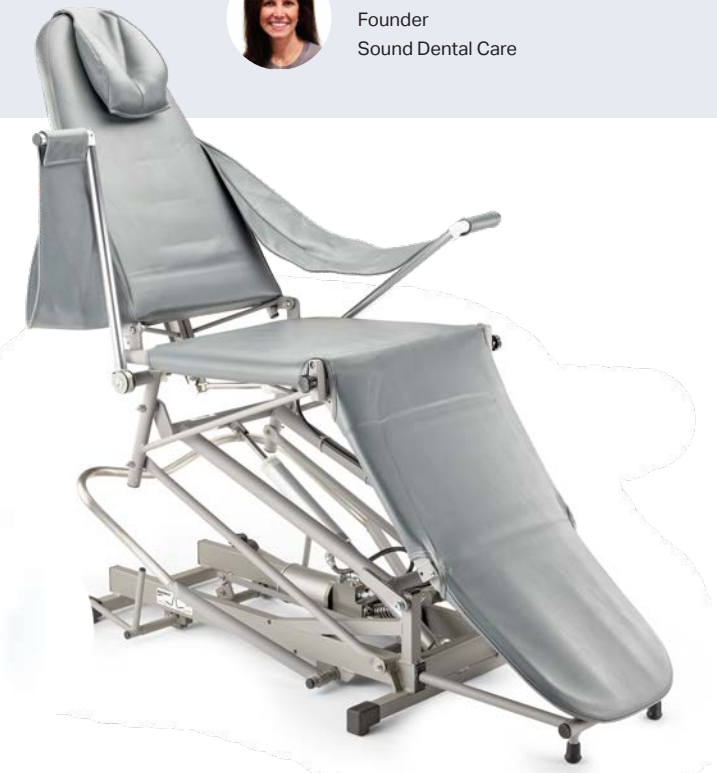

Hydraulic AseptiChair

ADC-02

- Hands-free height and back adjustment
- Adjustable with patient in chair
- 300 lb weight capacity
- No external electricity or air required
- Handle and wheels for easy mobility
- 69 lbs (31.3 kg)
- Hydraulic adjustable seat height: 11" - 33" (27.9 - 83.8 cm)
- Folded dimensions: 33" x 22" x 18" (83.8 x 55.9 x 45.7 cm)

The Hydraulic AseptiChair provides hands-free operation and is the only portable chair available with adjustability functions comparable to a standard dental chair. All chair lift functions are foot actuated hydraulic for hands-free operation; no external power is required. The AseptiChair provides comfortable patient seating for patients up to 300 lbs. Easy move wheels.

Hydraulic AseptiChair	ADC-02	\$6,990
Hydraulic Assembly Only	330629	\$1,800

Replacement Parts		
Replacement Dust Cover	410210	\$300



ALL ASEPTICO SEATING IS MADE IN THE USA. CALL TODAY! 425.408.9347

NOMAD™ PRO 2 Redefining Portable Dental X-ray Technology

ARU-07KP2

- Lightweight 5.5lbs (2.5kg)
- Simple touch controls on color LCD screen
- Removable and rechargeable battery packs
- One battery charge allows for more than 100 exposures
- Transport case included



The **NOMAD Pro 2 Portable Dental X-ray** offers cordless handheld operation, running on rechargeable battery packs.

Optimal for remote use, the NOMAD Pro 2 is contained in an easily transportable case with pull-out handle and wheels. Perfect for humanitarian & missionary work where portable dental x-ray is needed, or to conveniently move from station to station in a clinical setting.

Designed for operator safety and patient acceptance. Backscatter guard and internal shielding protect user.

Easy to operate and FDA approved.

Nomad Pro 2 Portable X-Ray System with Carrying Case ARU-07KP2 \$9,025

Replacement Parts

Premium Upgrade Kit

Includes Two Replacement Battery Handsets
One Replacement Charging Cradle
Power Cord

ARU-10134502 \$953

Basic Upgrade Kit

Includes Replacement Battery Handset
One Replacement Charging Cradle
Power Cord

ARU-10135062 \$589

ARU-10131887 \$519

ARU-10132534 \$270

Not Available In Canada



Portable LED Exam Light

ALU-40LED

- Adjustable LED with seven beam diameters
- Durable construction
- Includes padded carrying case
- Optional Amber UV filter for working with composites
- 12 lbs (5.4 kg)
- Padded carrying case: 45" x 9" x 5.5" (114.5 x 23 x 14 cm)

The **Portable LED Exam Light** provides high-intensity LED light with seven beam diameter settings. It features a compact lamp head with adjustable intensity, flexible gooseneck, telescoping tripod stand and padded carrying case. When a standard electric outlet is not available, an optional auxiliary DC input can be purchased that allows 12V and 24V battery operation in the field.

Portable LED Exam Light	ALU-40LED	\$1,530
Optional Amber UV Filter	461974	\$4
Replacement Padded Carrying Case	410084	\$160
DC Cable with Battery Clips	AA-40DC-BC	\$150
DC Cable with Car Lighter Adapter	AA-40DC-LP	\$170


Portable Tray Stand

ATC-03CF

The **Portable Tray Stand** is a sturdy and necessary addition to any mobile operator. Standard size stainless instrument tray removes easily for autoclaving.

It weighs less than 9 pounds and folds into included case when work is done. Made in USA.

Portable Tray Stand ATC-03CF \$960

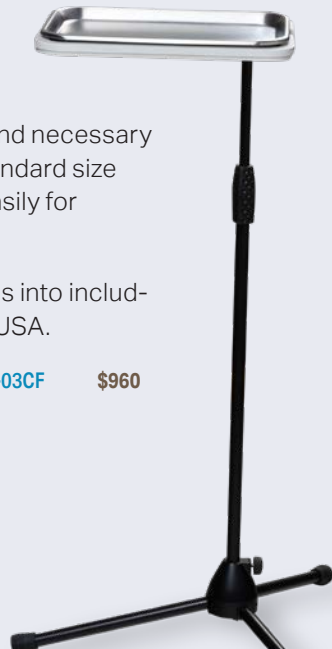

FRIENDLY AND PROFESSIONAL SERVICE. CALL TODAY! 425.408.9347



CALL TODAY! 425.408.9347

FOR FRIENDLY AND PROFESSIONAL CUSTOMER SERVICE.

Heavy Duty Brushless Lab Motor Torque Pro Delivers Power

AEU-16SS-120

- For heavy acrylic use
- Brushless 50k rpm micromotor
- Consistent torque even at low speeds
- Includes Variable Speed Foot Control



The Heavy Duty Brushless Lab Motor

delivers power, performance, durability, and portability. The brushless 50k rpm micromotor can produce up to 8.3 Ncm of torque and features easy to read digital display, preset memory storage, and hand/foot change button.

Heavy Duty Dental Lab Motor **AEU-16SS-120** **\$1,150**

Replacement Parts:

Motor/Cable Assy	AE-H180	\$540
Variable Speed Footswitch	AE-SP	\$75

Light Duty Lab Motor Powerful & Lightweight

AEU-10SS

- For light to medium acrylic use and low speed prophylaxis
- Versatile E-type straight 1:1 handpiece included
- 20k rpm micromotor
- Available in 110V or 220V
- Includes AE-7P On/Off Foot Switch



Includes AHP-101 Handpiece

The Light Duty Lab Motor

is one of the finest motor and handpiece combinations available to the professional. It is ideal for chairside adjustments, prophylaxis, or any light to medium use.

Light Duty Dental Lab Motor **AEU-10SS** **\$1,050**

Light Duty Dental Lab Motor **AEU-10SS-220** **\$1,050**

Options

Optional Variable Speed Foot Control **AE-14A** \$360

Replacement Parts:

1:1 Handpiece	AHP-101	\$380
Micromotor	AM-103B	\$220
On/Off Foot Switch	AE-7P	\$190

Battery Powered Portable Lab Motor Ideal for Mobile Hygiene

AEU-03M1B

- For light duty use and low speed prophylaxis
- Straight 1:1 e-type handpiece included
- Dual voltage 110V-240V



The Battery Powered Portable Lab Motor is a smart solution for lightweight lab benchwork and performing prophylaxis in the field! This small portable unit delivers 22k rpm and can perform 8 hours on a single battery charge.

Portable Lab Motor with Handpiece **AEU-03M1B** **\$800**

Portable Lab Motor **AEU-03M1B-NHP** **\$570**

Replacement Parts

1:1 Handpiece	AHP-101	\$380
Micromotor	AE-5M1-25	\$180
Foot Switch	AE-7M1	\$23

Portable Sterilizer

ASU-12-110 and ASU-12-230

- Compact and portable
- 18-minute sterilization cycle (24-minute total cycle time)
- Only 6 minutes to reach sterilization temp
- Temperature: 250°F (121°C)
- 11.5 lbs (5.2 kg)
- 9 liter capacity



The Portable Dental Autoclave is designed to sterilize solid, unwrapped instruments in a 121°C saturated steam process during 18-minute sterilization cycle. Use distilled or deionized water only.

110V Sterilizer **ASU-12-110** **\$1,399**

230V Sterilizer **ASU-12-230** **\$1,499**

Portable Saliva Ejector

AA-20

- Great for mobile dentistry or as a backup for sedation dentistry
- Powered by rechargeable battery pack, wall outlet or 12 VDC vehicle source
- Up to 60 minutes of continuous use on a single charge
- 800cc waste canister with float shut-off
- 6.3 lbs (2.9 kg)
- 9" x 7" x 8" (22.8 x 17.8 x 20 cm)



Great for mobile dentistry or as a backup for sedation dentistry, the Portable Saliva Ejector can be powered by a wall outlet, 12V source, or rechargeable battery pack. It features an adjustable 80 - 550 mm Hg vacuum range and 800cc waste canister with float shut-off. Includes padded carrying case.

Portable Saliva Ejector **AA-20** **\$590**

Replacement 800 ml Canister **AA-20CAN** **\$34**

Bacteria Filters **AA-20BF-2** **\$12**

Vacuum Pump Tube Assembly **330668** **\$49**

In addition to providing easy portable suction, the Portable Saliva Ejector fulfills the legal requirement for a back-up suction system in the event of power or suction failure.

Shop online at aseptico.com



CALL TODAY! 425.408.9347

FOR FRIENDLY AND PROFESSIONAL CUSTOMER SERVICE.

Customer Service 425.408.9988

Monday-Friday
6:00 a.m.-4:30 p.m. PST
Customer service representatives
are standing by to help.
repairs@aseptico.com

Sales 425.408.9347

Monday-Friday
6:00 a.m.-4:30 p.m. PST
Sales representatives
are ready to assist you.
orders@aseptico.com

Quality System

Aseptico's Quality Management System is certified to ISO 13485 as well as the Medical Devices Single Audit Program Recognized By The U.S. and Canada.

Return and Repair Policies

- Products that have been in clinical use cannot be returned for credit.
- Service and repair charges are billed on credit card. Customer is expected to accept charges up to \$750 without advance notification. Minimum service/repair charge is \$150.
- All returns and repairs require prior authorization. Contact Aseptico Customer Service for a Return Material Authorization number (RMA#).
- The RMA# must appear clearly on the exterior of the package being returned or it may be rejected and returned to the sender.
- Packages received without prior authorization are not the responsibility of Aseptico. Aseptico reserves the right to dispose of such packages and their contents 30 days after receipt with no obligation to the sender.
- The sender is responsible for proper packaging and for all costs for return shipment to Aseptico.
- New or unused products may be returned for full credit (less shipping charges) under a Return Material Authorization number within 30 days of date of original invoice.

Warranty

Aseptico warrants its products against defects in material or workmanship for a period of two (2) years, from date of original invoice. This warranty is not transferrable. Some handpieces are warranted for one year under the same conditions. Other handpieces and expendable components, such as air turbines and light bulbs, are covered by shorter warranty periods, or have no warranty. Aseptico's sole obligation under product warranty is (at its sole option and discretion) to repair or replace any defective component or product in part or whole. Aseptico shall be the sole arbiter of such action.

In the event of alleged defect under warranty, the purchaser is to notify Aseptico's Customer Service Department promptly. Customer Service will provide instructions, usually directing that the product be returned for service. Shipment to Aseptico and the cost thereof is always the responsibility of the purchaser.

Accidental misuse, inappropriate installation, or failure to perform directed maintenance voids the warranty.

Aseptico does not assume, under this warranty, any risks or liabilities arising from the clinical use of its products, whether or not such use involves coincidental utilization of products manufactured by others. Deliberately defacing, modifying, or removing the serial number voids the warranty.

Ordering Policies

- Aseptico accepts Visa, MasterCard, Discover, and American Express.
- Open account is available only if credit is approved in advance.
- Freight costs are added to all invoices.
- All orders are shipped FOB origin.
- Delinquent accounts are subject to 1.5% monthly interest, attorney's fees, collection and court costs when needed to satisfy debt.
- Prices subject to change without notice.
- Please visit: aseptico.com/warranty-policies for additional warranty and repair policy information.

Sales Tax

Aseptico uses a tax service (Avalara) to calculate and apply state sales tax to invoices for all states that Aseptico is registered to collect/remit sales tax (now including ID). We reserve the right to register in additional states at any time and begin collecting/remitting sales tax.

The amount of tax charged on your order will depend upon many factors including type of item purchased and destination of the shipment. Factors can change between the time you place an order and the time of invoicing, which could affect the calculation of sales taxes. The amount appearing on your order as sales tax may differ from the sales taxes ultimately charged.

Customer is responsible for any applicable state sales or use taxes not charged by Aseptico.

Sales Tax Exemptions


Unless Aseptico is in possession of a validated exemption certificate prior to shipment, the customer must pay the entire invoice, including sales tax. Certificates can be emailed to AR@aseptico.com. Check your Tax Exempt using our Avalara Cert Express Tool: Aseptico.com/tax-exempt

International Customers

Most Aseptico products are CE-marked and Health Canada licensed. However, customers should verify the status of any specific item(s) with their sales representative or local distributor prior to finalizing the purchase. Please contact your local Aseptico distributor for pricing and availability.

Prices shown in the catalog are in US dollars. If you are unsure of the Aseptico distributor in your country, for specific ordering and payment policies, or for general product inquiries, please contact our International Sales Department at int@aseptico.com.

Aseptico, Inc.
8333 216th St. S.E.
Woodinville, WA 98072

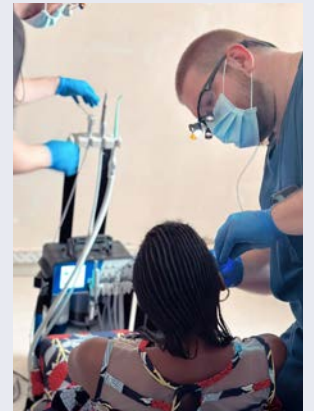

Giving Back Worldwide

Mission | Public Health | Schools | Nursing Homes | Military

Aseptico, Inc.
8333 216th Street SE,
Woodinville, WA 98072



Limited Time Package Offer*

15% OFF

Our Newest 1070 Series Implant Motor
& Mont Blanc Handpiece

*Exclusions Apply. Get 15% off the AEU-1070-70V and the AHP-85MB-X Handpiece. Cannot be combined with any other offer. Call 425-408-9347 for more details.

Welcome to the **Aseptico 2024 Global Dental Catalog**

Aseptico proudly manufactures products at its headquarters in Woodinville, WA. Aseptico® products are fully regulatory compliant to US and international manufacturing and safety standards, including FDA Good Manufacturing Practices, ISO 60601-1, and ISO 13485. These certifications, and our commitment to quality, allow us to export to over 125 countries.

425.408.9347

Aseptico.com



DENTAL MOTORS **WE POWER DENTISTRY**

MOBILE and PORTABLE DENTISTRY **MOVING TO MEET NEED**

GLOBAL DENTAL CATALOG **PROVEN PRODUCTS**

DENTAL EQUIPMENT OEM **EXCELLENCE IN ENGINEERING**