

Aseptico

2026

DENTAL CATALOG

NEW WIRELESS FOOT CONTROLS

Our Newest Innovation. Works with
Current Aseptico Motors (pg. 20)



PORTABLE DENTISTRY

Updated Taskforce Delivery System (pg. 31)

DENTAL MOTORS & HANDPIECES

Implant Motor: 1070 Series (pg. 4)



Welcome to the 2026 Aseptico Global Dental Catalog

As Aseptico moves into its 51st year, we can't help but be excited for what we have to deliver to the dental community at large in 2026.

Aseptico would like to introduce our latest innovation, the on/off and variable wireless foot controls. What we love about this is that it's not only an improvement to our modern systems, but is also available for many of our legacy motors as well, making this upgrade available to most of our customer base. Check out our spread for a list of compatible products!

It is not only the new that deserve to be talked about this year but also the improved, as one of our most popular portable dental systems, the ADU-12DCE Task Force, has an all-new design.

Additionally, we continue to make Aseptico a global enterprise as this past year we have achieved EU MDR certification for the AEU-350 GO Ultra Portable system. This means that our most advanced system is now available in many European markets.

You will find this year's catalog pages are filled with these new updates, plus our full line of motor systems and handpieces, innovative and trusted products such as HandiDam pre-framed dental dams, Penguin II Implant Stability Monitor, Spotter Smart Implant Detector, and much more.

We take pride that almost all Aseptico products are made in the USA and invite you to compare Aseptico products to any other for quality, performance, innovative features, reliability, and value.

We look forward to serving you in the coming year and beyond!

Glenn Kazen
President of Aseptico Inc.

Aseptico

DENTAL CATALOG



Dental Motors & Handpieces

WE POWER DENTISTRY

Recommended Systems	3
Oral Surgery and Implant Motors	3-8
1070 Series Implant Motor	4-5
Electrical Operative Motor	9
Endodontic Motors	10-11
Dental Handpieces	12-18
Tubing, Accessories & Replacement Parts	19-20



Innovations & Utility

PROVEN PRODUCTS

Penguin II Implant Stability Monitor	22
Spotter Smart Implant Detector	23
Powered Trolley Carts	26
HandiDam® Pre-Framed Dental Dam	27
Ortho Technique Chain & Grid	27
Dental Dam Clamps	27
Dental Dam Punches and Forceps	27
Dental Delivery Systems	28



Mobile & Portable Dentistry

MOVING TO MEET NEEDS

Recommended Systems	30
Portable Delivery Systems	31-37
Task Force Delivery System	31
GO Ultra Portable Delivery System	32-33
GO Cart Portable Supply Cart	34
Mobile Dental Carts	38-40
Portable Chairs & Stools	42-43
NOMAD Pro Portable X-ray	44
Portable LED Exam Light	45
Portable Tray Stand	45
Lab Motors	46
Portable Saliva Ejector	47

Customer Service & Ordering

VALUE & EXCELLENCE

Ordering Information, Policies and Contact Information	48
Repair and Product Support	48

Call 425.408.9347 or Shop Online at www.aseptico.com

Aseptico

DENTAL MOTORS AND HANDPIECES



AEU-1070
Our affordable implant motor streamlines procedures

AEU-7000
Award-winning, all-in-one premier surgical platform

AEU-6000
Our most popular and best value all-in-one platform

AEU-27
Excellence in rotary endodontics

	Recommended Use			
Oral Surgery		✓	✓	
Dental Implantology	✓	✓	✓	
Rotary Endodontics		✓	✓	✓
Restorative Dentistry		✓	✓	
	Motor and Handpiece Connection Type			
	Electric	Electric	Electric	Electric
	Special Features			
Wireless Foot Control Compatability	✓	✓	✓	✓
Multi-Function Foot Control	Yes or On/Off Option	Yes or On/Off Option	Yes or On/Off Option	
Handpiece Calibration	✓	Yes, Dynamometer	✓	✓
Handpiece Illumination		Yes, AEU-7000L-70V		
Dentsply Sirona rotary endodontic file library		✓		✓
User-customizable memory presets	✓	✓	✓	✓
Auto-stop-reverse for endodontics		✓	✓	✓
Upgradeable Software		✓	✓	✓
Autoclavable Tubing	✓	✓	✓	
	General Pricing*			
	\$3,300-\$3,500	\$5,230-\$5,750	\$3,730-\$3,950	\$2,300

Find The Right Dental Motor For You

Driven to advancing dental care and patient access to care. Call us today! 425-408-9347 or visit Aseptico.com

1070 Series Streamline & Simplify Implant Procedures

AEU-1070-70V

The Aseptico 1070 Series Implant Motor is designed to simplify dental implant procedures by providing all the features most desired by today's implant dentists with a straight forward design, smaller footprint and a more affordable price.



Watch a Video
Demonstration



Mont Blanc AHP-85MB-CX shown (not included). See details and other compatible handpieces pages 12-18.



Wireless Foot Controls Now Available

View Page 20 For More Information

Streamlined workflow

The 1070 Series Implant motor has the ability to place implants up to 70 Ncm with a 20:1 reduction handpiece and features a 1:1 setting for surgical and low-speed procedures.

User-friendly interface

The touchscreen display and intuitive user interface make it simple to adjust motor settings based on your preferences.

Hands-free control

The multi-function foot control enables easy adjustment of speed, direction, irrigation flow, torque, and presets, enhancing efficiency and infection control.

Precise and powerful

The brushless 40,000 RPM motor delivers the power needed for accurate and efficient implant procedures.

Versatility

With the 1:1 setting (and available 1:1 straight handpiece attachment), the motor is also suitable for surgical and low-speed procedures as well.

Customization options

The eight programmable presets allow you to customize the motor settings based on your preferred techniques.

Enhanced accuracy

The integrated handpiece calibration feature of the 1070 Series Implant Motor ensures speed and torque accuracy at the bur, providing consistent and reliable performance.

Compact and durable

The motor's compact footprint and water-resistant chassis design make it suitable for various clinical environments.

AEU-1070-70V

- Recommended for dental implants
- Up to 70 Ncm to streamline implant workflow
- Compatible with E-type 20:1 and 1:1 handpieces with external irrigation connections
- Advanced features and efficient adjustable torque desired by implant dentists
- Autoclavable 40krpm brushless micromotor
- Works for all implant types
- Adjustable irrigation flow rate of up to 140ml/min.
- Reusable, autoclavable irrigation tubing
- Easy cleaning and maintenance
- Console dimensions: 19.6cm W x 18.0cm D x 16.0cm H (7.72" W x 7.09" D x 6.30" H)
- Universal voltage 100-240V

NEW 1070 Motor with Wireless Multi-Function Foot Control	AEU-1070-70V-W	\$4,140
1070 Motor with Multi-Function Foot Control	AEU-1070-70V	\$3,700
1070 Motor with On-Off Foot Switch	AEU-1070	\$3,480

Options		
Micromotor w/ 10ft Cable	AE-230A-40EXT	\$1,370
Multi-Function Foot Control w/ 10ft Cable	AE-70V2-10	\$690
On/Off Foot Control w/ 10ft Cable	AE-7PM-10	\$240
Extended Irrigation Tubing Set	AE-24-10	\$120

Replacement Parts		
Micromotor w/ 6ft Cable	AE-230A-40	\$1,250
Multi-Function Foot Control	AE-70V2	\$620
On/Off Foot Control	AE-7PM	\$220
Irrigation Bag Hanger Rod	461541	\$48
Power Cord	840107	\$23
Handpiece Cradle	461561-01	\$11
Irrigation Tubing Set	AE-24	\$80
Irrigation Pump Tubing Section (qty: 10)	AE-24-PUMP	\$28
Tubing and Accessories see page 19-20		

7000 Series Motor All-in-One Surgical Platform

Oral Surgery, Implant Placement,
Endodontics & Restorative Dentistry

AEU-7000



AEU-7000L with Mont Blanc AHP-85MBFO-CX implant handpiece. Handpiece not included.
See details and other compatible handpieces pages 12-18.

The Aseptico 7000 Series Dental Motor

is one of the most powerful and easy-to-use implant/surgical motors available today. This all-in-one platform provides specialist-level performance for restorative, oral surgery, implant, and rotary endo applications.

All-in-One Platform

The 7000 Series has the ability to place implants up to 80 Ncm with a Mont Blanc® 20:1 handpiece.

For oral surgery, predictable and adjustable torque adds to practice efficiency.

For endodontics, the 7000 Series features auto-stop-reverse functionality to prevent file separation and comes preloaded with the complete Dentsply Sirona rotary file library plus the ability to customize your own presets.

For restorative procedures, the 7000 Series streamlines workflow; simply change out the handpiece and you have a restorative powerhouse.

Dynamometer Handpiece Calibration

The 7000 Series dynamometer calibration system gives dentists the ability to quickly and easily diagnose the health of a handpiece and verify that the handpiece is producing the correct torque for the application.

Upgradeable Software

Upgradeable software means a fantastic return on investment while ensuring that your future needs are also met.

Autoclavable Tubing

Reusable, autoclavable irrigation tubing means less waste, no more buying or keeping track of costly disposables.

Multi-Function Foot Control

The Multi-Function Foot Control replaces the need to adjust the unit manually during the procedure providing the operator variable speed control, forward-reverse, irrigation flow rate, torque adjustment, and memory functions.

Available handpiece ratios

1:1, 1:2, 1:3, 1:5, 8:1, 20:1

AEU-7000

- Recommended for oral surgery, dental implants, rotary endodontics & restorative dentistry
- All-in-one platform streamlines workflow
- Compatible with E-type 8:1 endo, 20:1 implant, 1:1, 1:2 surgical, 1:5, 1:3 high speed handpieces with external irrigation connections
- Additional handpiece ratios available via dynamometer calibration 4:1, 16:1, 32:1
- Advanced features and efficient adjustable torque desired by specialists
- Dynamometer provides exceptional handpiece testing and calibration
- Autoclavable 40k rpm brushless micromotor
- Works for all implant types
- Endodontic presets and complete Dentsply Sirona rotary endodontic file library onboard
- Upgradeable software
- Adjustable irrigation flow rate of up to 140 ml/min.
- Reusable, autoclavable irrigation tubing
- Easy cleaning and maintenance
- Optional: fiber-optic illumination, multifunction foot control, carrying case
- Console dimensions 9.98"W x 9.42"L x 5.1"H (25x24x13 cm)
- Universal voltage 100-240V

NEW 7000 Motor with Wireless Multi-Function Foot Control	AEU-7000E-70V-W	\$6,090
NEW 7000 Motor with LED Illumination and Wireless Multi-Function Foot Control	AEU-7000L-70V-W	\$6,390
7000 Motor with Multi-Function Foot Control	AEU-7000E-70V	\$5,650
7000 Motor with LED Illumination and Multi-Function Foot Control	AEU-7000L-70V	\$5,950
7000 Motor with On/Off Foot Control	AEU-7000E	\$5,430

Options

Carrying Case	410205-HNDL	\$950
Micromotor with 10ft cable	AE-230M-40EXT	\$1,420
Fiber Optic Micromotor with 10ft cable*	AE-230L-40EXT	\$1,470
On/Off Foot Control	AE-7PM	\$220
On/Off Foot Control with 10ft cable	AE-7PM-10	\$240
Multi-Function Foot Control with 10ft cable	AE-70V2-10	\$690
Extended Irrigation Tubing Set	AE-23-10	\$110

Replacement Parts

Fiber Optic Micromotor*	AE-230L-40	\$1,350
Micromotor with 6ft cable	AE-230M-40	\$1,300
Multi-Function Foot Control	AE-70V2	\$620
Calibration Bur	461558	\$16
Irrigation Bag Hanger Rod	461541	\$48
Power Cord	840079	\$23
Handpiece Cradle	461561-01	\$11
Irrigation Tubing Set	AE-23	\$80
Pump Section Tubing	AE-23-PUMP	\$28

More Tubing and Accessories see page 19-20

*Requires AEU-7000L Console



**WirelessFootControls
Now Available**

View Page 20 For More Information



"This is the best implant motor on the market. I've been using it for 8 years now and I've never had an issue. Very happy with this reliable equipment."

John Downey, DDS



Watch a Video
Demonstration



Blog: Aseptico® 7000 Series
Dental Motor is Named Top
25 Implant Product



PACKAGE DISCOUNTS AVAILABLE.

CALL TODAY! 425.408.9347

6000 Series Motor All-In-One Platform Our Best Value!

Osteotomy & Implant Placement,
Rotary Endo & Restorative

AEU-6000

- Recommended for implantology and rotary endodontics
- Wireless foot controls now available
- All-in-one platform streamlines workflow
- Compatible with E-type 8:1 endo, 20:1 implant, 1:1, 1:2 surgical, 1:5, 1:3 high speed handpieces with external irrigation connections
- Integrated Handpiece Calibration
- Autoclavable 40k rpm brushless micromotor
- Works for all implant types
- Endodontic presets
- Upgradeable software
- Adjustable irrigation flow rate of up to 140 ml/min.
- Reusable, autoclavable irrigation tubing.
- Easy cleaning and maintenance
- Small footprint:
8.6 W x 9.2" L x 4.8" H (22 x 23 x 12 cm)
- Universal voltage: 100-240V



AEU-6000 shown with On/Off Foot Pedal and Mont Blanc Handpiece (Not Included)

"I have really enjoyed my Aseptico AEU-6000. It meets all my surgery needs from placing implants to removing full bony impacted wisdom teeth. I have found that this unit has more features and better torque than my old W&H."

Heath Snell, DDS
General Dentist

The AEU-6000 is our best-selling all-in-one motor because of its reliability and affordability. This versatile machine has the ability to place implants up to 50 Ncm with a Mont Blanc 20:1 handpiece. It performs rotary endo with auto-stop-reverse, and it can also be used as a backup for your restorative and oral surgery systems. No other motor on the market can give you these benefits at such an affordable price!

With upgradeable software, the AEU-6000 unit will be a workhorse for many years to come. Autoclavable irrigation tubing means you don't have to buy a new set for each procedure. Five preset buttons can be programmed separately for both implants and endodontics, reducing your procedure time dramatically.

Available handpiece ratios
1:1, 1:2, 1:3, 1:5, 8:1, 20:1

NEW Implant Motor w/ Wireless Multi-Function Foot Control	AEU-6000-70V-W	\$4,590
NEW Implant Motor w/ Wireless On/Off Foot Switch	AEU-6000-W	\$4,340
Implant Motor w/ Multi-Function Foot Control	AEU-6000-70V	\$4,150
Implant Motor w/ On/Off Foot Switch	AEU-6000	\$3,930

Options		
Carrying Case	410205-HNDL	\$950
Micromotor with 10ft cable	AE-230M-40EXT	\$1,420
On/Off Foot Control	AE-7PM	\$220
On/Off Foot Control with 10ft cable	AE-7PM-10	\$240
Multi-Function Foot Control with 10ft cable	AE-70V2-10	\$690
Extended Irrigation Tubing Set	AE-23-10	\$110

Replacement Parts		
Micromotor with 6ft cable	AE-230M-40	\$1,300
Multi-Function Foot Control	AE-70V2	\$620
Irrigation Bag Hanger Rod	461541	\$48
Power Cord	840079	\$23
Handpiece Cradle	461561-01	\$11
Irrigation Tubing Set	AE-23	\$80
Pump Section Tubing	AE-23-PUMP	\$28

Tubing and Accessories see page 19-20

Electrical Operative Motor Top-Rated Motor Converts Air to Electric

AEU-5000



Recommended handpiece 1:5 Restorative with fiber-optic illumination (AHP-72MK-FO)

Available handpiece ratios 1:5, 1:2, 1:1, 5:1, 8:1

The versatile and cost-effective operative motor is designed to convert a standard air dental delivery system into a powerful electric handpiece system for numerous dental applications.

The AEU-5000 features endo functionality, fiber optic intensity adjustment, a compact and easy-to-clean console, a large backlit display, intuitive user interface, and a 40,000 rpm brushless autoclavable micromotor. Innovative LED technology produces some of the brightest illumination you'll find.

- Converts air to electric
- Powerful, high torque electric motor
- Attach to your air-driven high speed connection and use your rheostat
- Electrical option for efficient and quiet restorative procedures
- Up to 200,000 rpm with high speed handpiece
- Intuitive design and operation
- Includes endo mode with auto-stop-reverse torque setting
- Universal voltage 100-240V

Electrical Operative Motor	AEU-5000	\$1,990
Adjustable Mounting Bracket	330569	\$110
Replacement Micromotor	AE-240-40	\$1,750
1:5 Restorative Handpiece with Fiber-Optics	AHP-72MK-FO	\$1,070

Recommended Handpieces see pages 12-18



Optional AMC-25 Integration
Read more about the benefits of our most advanced dental cart on page 38.

AMC-25 Product Shown Not Included, Sold Separately

Endo DTC® Rotary Motor America's Most Popular

AEU-25

- Most sold, most trusted, endo motor in America
- Dentsply Sirona file library included
- Auto-stop-reverse functionality



Foot control included



MicroNiTi Handpiece

We recommend the AHP-88MNP MicroNiTi endodontic handpiece. See page 15.

"As a busy endodontist, I want equipment that is reliable, but more importantly, from a company that has great customer service. Aseptico delivers on both fronts."

Chad Burgess, DDS
Endodontist

The Endo DTC Rotary Endo Motor is America's most popular endo motor because it is a proven reliable performer. Many dentists did their endodontic training on this system and they still prefer it today. There is so much packed into this small footprint at such an affordable price: advanced digital torque control technology with auto-stop-reverse helps prevent file separation, five programmable memory buttons, onboard file library, integrated calibration system, autoclavable 30k rpm brushless micromotor, and upgradeable software. Pair with a AHP-88MNP handpiece to make the perfect rotary endo system.

Available handpiece ratios
1:5, 1:1, 4:1, 8:1, 16:1, 64:1

Endo DTC® **AEU-25** \$2,250

Replacement Parts

Micromotor	AE-4B-30	\$850
On/Off Foot Control	AE-7P	\$220
Handpiece Holder	AE-5FG	\$26
Motor Cap (For AE-4B-10/30)	461084-08	\$10
Motor Plug (For AE-4B-10/30)	460183-08	\$11

Recommended handpieces see page 15



Electronic Endodontic System Our Most Advanced Endo Motor

AEU-27A

- Lightweight motor and cable
- Upgradeable software
- Weighted omnidirectional foot pedal
- Includes all features of America's most popular endo motor, the Endo DTC®



Don't forget the HandiDam®

Created by Aseptico, it's the original pre-framed rubber dam trusted worldwide. Easier and more efficient — it's the only rubber dam that allows access to suction or take x-rays while isolation and clamps remain in place. See page 27.



The Electronic Endodontic System features digital torque control with auto-stop-reverse to prevent file separation and automatic calibration to test both the motor and handpiece for accurate performance and consistent operation. Easy to operate, with an auto-clavable brushless motor, this motor also features onboard file library, five programmable memory locations, and preset options that allow single file, file series, or custom file series. Upgradeable software assures return on investment. Includes weighted omnidirectional foot pedal.

Available handpiece ratios
1:5, 1:1, 8:1, 16:1

Electronic Endodontic System **AEU-27A** \$2,450

Options

NEW Wireless Foot Switch **AE-7PMW-4** \$620

Replacement Parts

Micromotor	AE-250-30	\$730
On/Off Foot Switch	AE-7PM	\$220
Motor Cap (For AE-4B-10/30)	461084-08	\$10
Motor Plug (For AE-4B-10/30)	460183-08	\$11

MicroNiTi endodontic handpiece recommended see page 15

"The AEU-27A is the most versatile, safe and economical motor available anywhere in the world, at any price. Upgradeable, simple to use, and ergonomic, it sets the standard for electric endodontic motors for many years to come."

Richard Mounce, DDS

CALL TODAY! 425.408.9347 or shop online at aseptico.com

Mont Blanc® World Class Implant Handpieces

The Mont Blanc Implant Handpiece is crafted from stainless steel and titanium, making it lightweight and durable. Precision ball bearing construction reduces friction to prevent overheating during surgical procedures. It features smooth surfaces that improve infection control and ensure easy cleaning, with the ability to remove the head for cleaning without tools.

Electric (E-type) | 80 Ncm max torque | 20:1 reducer | 15-2,000 rpm (40k rpm micromotor)



Economy Best Value Implant Handpiece

With an upper torque limit of 70 Ncm, the Economy Dental Implant Handpiece is ideally suited as a primary or backup handpiece for the AEU-1070 and AEU-6000 platforms. This implant handpiece offers exceptional value and is excellent for drilling osteotomies, tapping and threading.

Electric (E-type) | 70 Ncm max torque | 20:1 reducer | 15-2,000 rpm



AseptiSpray™ Ensuring Handpiece Longevity



Enhance Your Handpiece Performance and Longevity with AseptiSpray™ Advanced Cleaner and Lubricant:

- Provides complete internal care for all types of air-driven (pneumatic) and E-type (electric) handpieces.
- Made by Aseptico to ensure the highest standards for modern dental equipment.
- Features a CFC-free and ozone-safe formula for responsible practice.
- Proper lubrication enhances performance and extends the life of your handpieces.

Learn More On Page 19.

Watch a Video
Demonstration



ASEPTICO SELECT HANDPIECES ARE EXCEPTIONAL PRODUCTS,
RIGOROUSLY QUALIFIED AND PERFORMANCE TESTED.

See all handpiece accessories, including AseptiSpray™ handpiece cleaner and lubricant, irrigation clips, cannulas, and more on pages 19-20.

20° Surgical

AHP-77W

- Electric (E-type)
- 20° bend for optimum comfort
- 1:2 increaser
 - 600 to 80,000 rpm (w/ 40k rpm micromotor)
 - Requires 70 mm bur length



AHP-77W 20° Electric High Speed Handpiece is the oral surgeon's number one choice for third molar extractions. The angled design prevents fatigue by allowing it to rest on the back of the hand. A scratch-resistant surface and reinforced nose bearings guarantee durability, and the sliding lip seal protects the interior from penetration of dirt particles.

20° Electric Handpiece	AHP-77W	\$2,627
W&H Cleaner and Lubricant	AHP-10 & AHP-10N	\$92
Irrigation Clip	AHP-77W-CLP	\$86

45° Surgical

AHP-79W

- Electric (E-type)
- 1:2.7 increaser
- 810 to 108,000 rpm (w/ 40k rpm micromotor)



W&H 45° High Speed Surgical Handpiece is ideal for precise tooth separations and apical resections. The 45° head allows for optimum access to the posterior of the mouth. Available in non-fiber optic and self-generating fiber-optic options, this handpiece comes complete with a triple-spray head that guarantees sufficient cooling of the surgical bur and site.

45° Electric Handpiece w/LED	AHP-79W-FO	\$2,724
45° Electric Handpiece	AHP-79W	\$2,261
W&H Cleaner and Lubricant	AHP-10 & AHP-10N	\$92

1:2 Straight Surgical

AHP-77DS

- Electric (E-type)
- 1:2 increaser
- Up to 80Krpm



AHP-77DS Electric High Speed Handpiece is lightweight and compact. The traditional straight design allows for exceptional posterior access. With its 1:2 increaser and rotation speeds up to 80,000 rpm, this handpiece is a great option for bone contouring and sectioning teeth.

1:2	AHP-77DS	\$2,150
Replacement Irrigation Clip	AHP-77D-NOZ	\$34
AseptiSpray™ Cleaner and Lubricant	AHP-09	\$46

Microhead Endodontic

AHP-88MNP

- Electric (E-type)
- Push button
- Speed range (w/ 30k rpm micromotor)
- 8:1 — 40 to 3,600 rpm



Anthogyr MicroNiTi Microhead Endodontic Handpiece has a push button head and provides the same unrivaled accuracy as other MicroNiTi series handpieces, whether it is being used with a rotary or reciprocating motor. The micro-head size allows for exceptional access and line of sight during procedures. Made for nickel-titanium endodontic files, the MicroNiTi Microhead is compatible with all ISO rotary endo files.

8:1	AHP-88MNP	\$830
AseptiSpray™ Cleaner and Lubricant	AHP-09	\$46

1:3 Restorative

AHP-65TI-FO and AHP-65TI

- Electric (E-type)
- External irrigation
- 1:3 increaser
- Optional fiber optics
- Max 120,000 rpm (w/ 40k rpm micromotor)
- Perfect for high speed dentistry with your implant/surgical motor



With an external irrigation port and 1:3 increaser gearing, the AHP-65TI Restorative Electric High Speed Handpiece makes any surgical motor with irrigation into a restorative powerhouse. Ergonomic and lightweight design allows for optimum comfort during your procedure.

1:3 High Speed w/ Fiber Optics	AHP-65TI-FO	\$1,800
1:3 High Speed	AHP-65TI	\$1,300
AseptiSpray™ Cleaner and Lubricant	AHP-09	\$46

1:5 Restorative

AHP-72MK-FO

- * Electric (E-type)
- Internal irrigation
- 1:5 increaser
- Fiber optics
- Max 200,000 rpm (w/ 40k rpm micromotor)
- Perfect for high speed dentistry with your electric operative motor



Provides powerful, quiet, vibration-free operation with class-leading precision, the 1:5 Restorative Handpiece features Long Lasting Gear that are engineered to excel at every angle, and ceramic bearings that stand up to heavy use. Fiber optics optimize visibility with 100% glare-free illumination at a color temperature of 5500 K.

1:5 High Speed	AHP-72MK-FO	\$1,070
AseptiSpray™ Cleaner and Lubricant	AHP-09	\$46

4-Hole Restorative

AHP-96MK

- Air-driven
- Internal irrigation
- Non-optic
- Max 320k rpm, 24 watts
- Max volume 62 dB
- Recommended for high speed dentistry



Simple to use technology, developed for what you do every day. The 4-hole High Speed Handpiece features a clean, elegant shape and a virtually scratch proof smooth chrome surface. Manufactured using proprietary laser balancing technology, this handpiece delivers 24 watts of torque and features a premium cartridge for enduring precision and prolonged life.

MK Dent 4-hole High Speed Restorative Handpiece	AHP-96MK	\$540
AseptiSpray™ Cleaner and Lubricant	AHP-09	\$46

5-Hole Restorative

AHP-97MK

- Air-driven
- Internal irrigation
- Fiber optic glare-free, color accurate illumination
- Max 320k rpm, 24 watts
- Max volume 62 dB
- Perfect for high speed dentistry



It feels great in your hand and sports a clean modern shape with a chrome surface. Fiber optics deliver 100% glare-free, color accurate illumination. The 5-hole High Speed Restorative Handpiece is equipped with the AeroPerformance Blade premium cartridge. Ceramic bearings stand up to heavy use and autoclaving.

MK Dent 5-hole High Speed Restorative Handpiece	AHP-97MK	\$570
AseptiSpray™ Cleaner and Lubricant	AHP-09	\$46

KaVo-style Restorative

AHP-82MK

- Air-driven
- Internal irrigation
- Fiber optic with LED color accurate illumination
- Max 350k rpm, 25 watts
- Max volume 55 dB



This KaVo-style swivel restorative handpiece is ideal for high-speed dentistry in the office or with a mobile pneumatic delivery system. It delivers more power, quieter operation, and excellent control. The MK-Dent 6-pin High Speed Restorative Handpiece is built for professionals who demand top performance. Reverse-flow protection helps reduce contaminants entering the head, and fiber optics provide clear, color-accurate lighting. Finished with a scratch-resistant titanium surface.

Requires Coupler (AHP-82MK-QD or AHP-82EGEN-QD).

MK Dent 6-pin High Speed Restorative Handpiece	AHP-82MK	\$1,130
MK Dent 6-pin Coupler	AHP-82MK-QD	\$520
4 Hole Fiber Optic Coupler	AHP-82EGEN-QD	\$295
AseptiSpray™ Cleaner and Lubricant	AHP-09	\$46

6-Pin Fiber Optic Coupler

AHP-82MK-QD

- Air-driven
- Works with LED bulbs and glass rod fiber optics
- Color accurate



This 6-Pin Fiber Optic Coupler comes with a non-retraction valve and water volume regulation. Works with both LED and glass rod fiber optics to provide 100% glare free and color accurate light.

6-Pin Fiber Optic Coupler [AHP-82MK-QD](#) \$520

4-Hole Autogenerating LED Coupler

AHP-82EGEN-QD

- Air-driven
- Compatible with AHP-82MK



This innovative 4-hole coupler adds LED functionality to standard 4-hole handpiece connections without the need to connect to an electricity source. Recommended for Aseptico's air-driven economy air, portable, and mobile delivery units.

4 Hole Fiber Optic Coupler [AHP-82EGEN-QD](#) \$295



AHP-82MK Handpiece shown with AHP-82-MK-QD Coupler attached.

Low Speed Straight

AHP-101

- Electric (E-type)
- Attaches to E-type air motor
- Also pairs with AHP-50MK
- 1:1 Straight Handpiece
- Post for Doriot-style prophy angle attachment
- Recommended for use with Aseptico portable systems and lab motors



The Low Speed Straight Handpiece is made of lightweight aluminum and is designed to work with Aseptico® air-driven portable dental systems. The AHP-50 Low Speed Air Motor, which allows air turbine to accept E-type handpieces, is required to connect this handpiece to Aseptico portable dental systems.

Low Speed Straight Handpiece **AHP-101** \$380

Low Speed Contra Angle

AHP-63MK

- Electric (E-type)
- Attaches to E-type air motor
- Also pairs with AHP-50MK
- Swing Latch
- 1:1 Contra Angle
- Recommended for low speed work on portable equipment or in an office operator



Low Speed Contra Angle Handpiece is a D-latch style low speed handpiece that can be attached to any E-type air motor. Made by MK Dent

MK Dent 1:1 Swing Latch Contra Angle **AHP-63MK** \$290

AseptiSpray™ Cleaner and Lubricant **AHP-09** \$46

4-Hole to E-Type Air Motor

AHP-50MK

- Air Driven
- Low Speed Air Motor (4-hole)
- Adapts an air system to accept an E-type handpiece



The Low Speed Air Motor accepts E-type handpieces. Recommended for use with Aseptico® air-driven portable dental systems or other in-office applications where you need to adapt an air system to an E-type handpiece.

Low Speed Air Motor **AHP-50MK** \$410



AseptiSpray™ Handpiece Cleaner and Lubricant

AseptiSpray™ is a must for proper cleaning and maintenance of handpieces. It keeps handpieces working efficiently and ensures long life and durability.

AseptiSpray™ w/ Nozzle (8.8 oz can)	AHP-09	\$46
Nozzle only	AHP-09N	\$15

W&H f1 Service Oil Handpiece Cleaner and Lubricant

Synthetic oil of maximum purity, which has been specially made for maintenance of W&H handpieces. Required by W&H warranty policy.

Handpiece Service Oil (Spray Cap NOT included)	AHP-10	\$77
Service Oil Spray Cap	AHP-10N	\$15

Blog: Advantages of Autoclavable Tubing for Oral Surgery Systems



Mont Blanc®/Economy Handpiece Accessories

Depth Stop Attachment	AHP-85MB-STP	\$130
Black Irrigation Clip	AHP-85MB-CLP	\$60
Cannula	AHP-85MB-CAN	\$33
Cannula Clip Set (Includes Small Clip, Large Clip, Cannula and Tubing)	AHP-07	\$24
Cannula and Tubing	AHP-04	\$19
Small Clip	AHP-05	\$8
Large Clip	AHP-06	\$8
Y Connector	AE-2	\$9



Autoclavable Black Irrigation Clip for Mont Blanc® Handpiece **AHP-85MB-CLP**

Tubing & Accessories

Autoclavable surgical tubing saves money and prevents waste. Disposable tubing sets cost about \$12 per set compared to less than \$1 per use for autoclavable tubing that lasts for approximately 100 uses. Autoclavable irrigation clips are money-savers too.

Irrigation Tubing Set for IV Bags (For AEU-6000 and 7000 series)	AE-23	\$80
Extended Irrigation Tubing Set	AE-23-10	\$110
Irrigation Pump Tubing Section (qty:10) (For AEU-6000 and 7000 series)	AE-23-PUMP	\$28
Irrigation Clips (qty: 10)	AHP-03	\$22
Irrigation Tubing Set for IV Bags (For AEU-1000 and 1070 series)	AE-24	\$80
Extended Irrigation Tubing Set	AE-24-10	\$120
Irrigation Pump Tubing Section (qty:10) (For AEU-1000 and 1070 series)	AE-24-PUMP	\$28
Motor Cap (For AE-4B-10/30)	461084-08	\$10
Motor Plug (For AE-4B-10/30)	460183-08	\$11

CALL TODAY! 425.408.9347 or shop online at aseptico.com

FOOT CONTROLS

WIRELESS OPTION NOW AVAILABLE



Multi-Function Foot Control

AE-70V2

The Multi-Function Foot Control replaces the need to adjust the unit manually during procedures, providing the operator with variable speed control or on/off pedal, forward-reverse, torque adjustment, irrigation flow rate, and memory functions for implants and endodontics. Made in USA.


NEW Wireless Multi-Function Foot Control	AE-70V2W-4	\$1,060
Multi-Function Foot Control with 6ft cable	AE-70V2	\$620
Multi-Function Foot Control with 10ft cable	AE-70V2-10	\$690

On/Off Foot Switch

AE-7PM

Wired and Wireless On/Off Foot Switch works on all versions of AEU-26, AEU-27, AEU-28, AEU-1000 series, AEU-6000 series, AEU-7000 series, and AEU-350/S. Made in USA.

NEW Wireless Foot Switch	AE-7PMW-4	\$620
Foot Switch with 6' cable	AE-7PM	\$220
Foot Switch with 10' cable	AE-7PM-10	\$240

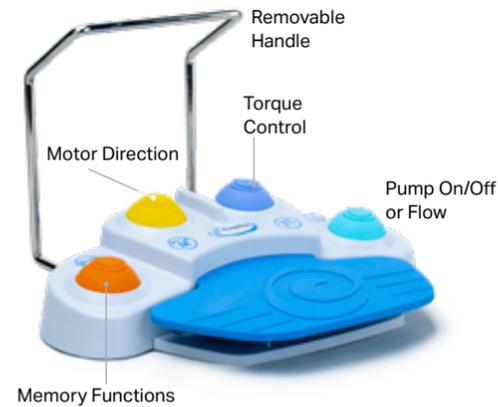

Replacement Micromotors

AEU-7000L micromotor and 6' fiber optic cable	AE-230L-40	\$1,350
10 ft extended cable option	AE-230L-40EXT	\$1,470
6000 Series and 7000 Series Implant Motors micromotor and 6' cable	AE-230M-40	\$1,300
10 ft extended cable option	AE-230M-40EXT	\$1,420
AEU-5000 Electric Operative Motor micromotor	AE-240-40	\$1,750
Micromotor and Cable for:	AE-4B-30	\$850
AEU-25 Endo DTC®		
AEU-707AV2 Implant/Surgical Motor		
AEU-17BV2 SurgiMotor II		

Motor Cap for AE-4B-10/30	461084-08	\$10
AEU-10SS Light Duty Lab Motor replacement micromotor	AM-103	\$230
AEU-1000/1070 Micromotor	AE-230A-40	\$1,250
AEU-1000/1070 Micromotor with 10 ft Cable	AE-230A-40EXT	\$1,370
AEU-27A Advanced Endodontic Motor micromotor	AE-250-30	\$730

Aseptico

INNOVATIONS
AND UTILITY





Newest Generation Penguin II Implant Stability Monitor

RFA-55250-MULTI

New and Better! Aseptico is proud to announce the US arrival of the next generation RFA implant stability monitor – the Penguin II.

- Evolved design brings increased comfort, surer grip
- Relocated, easier-to-read LED display screen
- Battery can now be replaced, extending product life
- New vertical wireless charging stand saves space

Spotter Smart Implant Detector Find Buried Implants Faster and Preserve Soft Tissue

SPT-5I11001-KIT

- Painlessly and accurately locate buried implants quickly and easily
- Non-invasive, no exploratory flap needed
- Minimize the recovery time for the patient
- Two-year warranty
- Works with any titanium implant

Spotter Smart Implant Detector locates dental implants covered by gum tissue quickly and painlessly — no exploratory flap needed! Just move the wand over the site and the Spotter will guide you to the exact center of the implant.

Instead of opening a flap to locate the implant, the Spotter finds the exact center quickly, allowing you to punch and reveal only the implant, thereby helping to preserve soft tissue and minimize recovery time for the patient.

The Spotter is especially helpful in the case of multiple implants when implant sites may be hard to visualize. The Spotter Smart Implant Detector works with any type of titanium implant.

It is easy to use, requiring only a few minutes to learn how to operate. The unit operates on one AA battery and its detachable sensing head is autoclavable.

Spotter SPT-5I11001-KIT \$700
Sensors 5/PKG SPT-5I11101 \$130



"With Spotter, we don't have to raise a flap, we can go in a flapless procedure, just punching the soft tissue to get to the center of the implant."

Ziv Mazor, DMD

spotter
Smart Implant Detector

Implant Stability Monitors

- Provides accurate Implant Stability Quotient measurement for greater confidence determining when to load dental implants
- Cordless, lightweight, and simple to use
- Half the price of competing ISQ devices
- Two-year warranty, Made in Sweden

Penguin Implant Stability Monitors

Now every dentist can access affordable Resonance Frequency Analysis (RFA) technology with the Penguin RFA and Penguin II for assessing osseointegration of dental implants. These devices measure local bone density, translating it into an implant stability quotient (ISQ) value between 1 and 99. Higher ISQ values indicate greater stability. Especially useful for all-on-4 or all-on-X cases, the Penguin monitors provide precise, objective insights into osseointegration, ensuring predictable outcomes.

MulTipeg™ Sensing Devices

Both Penguin RFA and Penguin II utilize the same MulTipeg™

- For all major implant systems
- Made of durable, tissue friendly titanium
- Autoclavable for 20 uses
- Laser marked with type number
- ISQ Standard Calibrated™
- Cost effective

Penguin II	RFA-55250-MULTI	\$2,690
Penguin RFA	RFA-55002-MULTI	\$2,440
Penguin MulTipeg	RFA-550**	\$42
MulTipeg Driver	RFA-55003	\$27
Sterile Covers (20 Pack)	RFA-55105	\$55
Penguin ISQ Tester	RFA-55217	\$16
Penguin II Replacement Battery	RFA-55291	\$210

**Part number varies by implant platform.

Watch a Video Demonstration



View MulTipeg List



CALL TODAY! 425.408.9347 or shop online at aseptico.com

penguin II WANT TO REDUCE IMPLANT FAILURE?

THE NEXT GENERATION OF
RFA TECHNOLOGY: **PENGUIN II**

Implant Stability Monitor now with replaceable
batteries and a charging stand

MEASURE IMPLANT
STABILITY IN SECONDS

For more info, call 425-408-9347
or visit Aseptico.com/Penguin-II



"I truly feel that the use of RFA
technology is imperative in
order to provide timely and
predictable implant treatment
to my patients."

Charles D. Schlesinger
DDS, FICOI



NEW
AEU-350S

Featuring
**ULTRASONIC
SCALER**

Aseptico
Trusted Worldwide

MOBILE & PORTABLE SOLUTIONS Dental clinicians, public
health outreach programs, humanitarian dental groups, and
the US military rely on Aseptico® mobile and portable dental
equipment to deliver dental care around the world. **Made in USA.**

Dental Care Anywhere

Aseptico® Mobile & Portable Solutions

- Mobile Dental Carts
- Portable Delivery Systems
- Patient & Clinician Seating
- Portable X-ray
- Portable Exam Lights



Portable Seating Up to Any Task!

FRIENDLY & PROFESSIONAL CUSTOMER SERVICE CALL TODAY! 425.408.9347
or shop online at aseptico.com



Aseptico Trusted Worldwide



Powered Trolleys

ATC-12V2 and ATC-15

- Keeps equipment organized and solves operatory space concerns
- Storage bin for foot pedals or other items
- Dual outlet with coiled AC power cord
- Adjustable height to easily store under a counter: 26.4" to 37.4" (67-95 cm)
- Four shelves w/ 3 removable stainless steel trays

Aseptic Powered Trolleys give dental teams the ability to easily roll dental motors, scalers, lasers, surgical systems or other electrical equipment from operatory to operatory. Available in large and standard sizes, each Powered Trolley is adjustable in height and features four shelves of storage with removable stainless-steel trays that fit securely into each shelf. High-quality, sturdy trolleys move and turn with ease and fit into almost any small space.

Standard Powered Trolley **ATC-12V2** \$845

Cart size: 26.4-37.4 x 15.75 x 16.5 (67-95 x 40 x 42 cm)
Tray size: 11.3" x 9" (28.8 x 22.8 cm)

Large Powered Trolley **ATC-15** \$1,000

Cart size: 26.5-38 x 19.5 x 19.75 (68-97 x 50 x 50 cm)
Large Tray: 15.4" x 11" (39.1 x 28 cm)



Aseptic Powered Trolleys feature two heavy duty electrical outlets that are certified to safety standards and a sturdy storage bin designed to hold foot pedal or other items.



HandiDam® pre-framed dental dam



- Pre-framed, flexible design facilitates access to the oral cavity for suction or X-rays
- Superior elasticity and tear-resistance
- Saves time—eliminates the need to remove and replace traditional dam during the procedure
- Ergonomic and comfortable for greater patient acceptance
- Meets FDA Powder-Free standard
- Extremely low protein content reduces patient irritation (<50 micrograms)
- Available in Latex or Latex-Free



HandiDam® Latex (Blue)

Qty 20, **I-70BM** \$45
Qty 100, **I-70BM-100** \$170

HandiDam® Latex-Free (Mauve)

Qty 20, **I-70LF** \$70
Qty 100, **I-70LF-100** \$300


TECHNIQUE CHAIN & GRID



Technique Chain & Grid easily attaches to a pre-emerged tooth, so that an orthodontic device may be used to draw the tooth into proper position. Aseptico is the originator of this remarkable traction system that professionals use worldwide.

- 10 sets per box
- 14k gold chain
- 3 grid shapes available

Cuspid (qty: 10) **AMA-09A** \$475
Bicuspid (qty: 10) **AMA-09A2** \$475
Incisor (qty: 10) **AMA-09A3** \$475

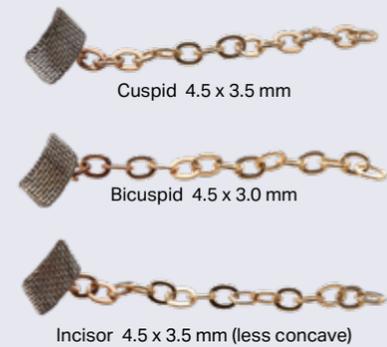

Dental Dam Clamps

Our five most popular clamps. See full selection at aseptico.com
Dental Dam Clamp \$23 each



K-211
Universal labial cavities on anteriors



K-5
For lower molars



K-9
Most popular anterior clamps



K-90N
Most popular anterior clamps



K-8A
For lower molars, with strong festooned beaks





Automated Air Delivery System

ADU-07 and ADU-07A

- Choose 2 or 3 handpiece system
- Optional fiber-optics
- Includes 3-way air/water syringe
- Includes foot control
- Wall or cabinet mount

The Automated Air Dental Delivery System is a pneumatic handpiece control system that gives you everything you need at a price you can afford. Highly customizable, select from a 2-handpiece system with manual control, or 3-handpiece system with automatic selection and individual air and water adjustments. Both models include air pressure gauge, 3-way air/water syringe, foot control, and water flow adjustments. Flat top design is ideal for mounting into cabinets. Requires connection to air and water sources.

2-Handpiece/Manual Control	ADU-07	\$1,200
3-Handpiece/Automatic	ADU-07A	\$1,250

Options

5-hole Fiber Optic Upgrade	AA-19LED-04	\$720
6-pin Fiber Optic Upgrade	AA-19A-06	\$720
4 Hole Fiber Optic Coupler	AHP-82EGEN-QD	\$295

Replacement Parts

Disc Foot Control w/ Water Toggle	AA-43W	\$110
Air/Water Syringe	TA-90D	\$85
Autoclavable Tips for Air/Water Syringe	TA-1	\$9



Folding Arm Wall Mount
Folding wall mount for ADU-07, ADU-07A and ADU-04, Dental Delivery Systems
AA-05A \$450



Manual Air Delivery System

ADU-04

- 2-handpiece control system
- Ideal for arm or wall mounting
- 3-way air/water syringe
- Requires connection to air
- Optional independent water

The Manual Air Delivery System is a pneumatic handpiece control system for two handpieces. It includes a 3-way air/water syringe and disk foot control. Requires connection to air and water sources or use optional independent water.

Manual Air Control System	ADU-04	\$680
---------------------------	------------------------	-------

Options

Clean Water System with 1L Bottle (Shown)	NWS-7I	\$130
4 Hole Fiber Optic Coupler	AHP-82EGEN-QD	\$295

Replacement Parts

Disc Foot Control	AA-42	\$65
Air/Water Syringe	TA-90D	\$85
1L Water Bottle (for NWS-7I)	730471	\$13
Autoclavable Tips for Air/Water Syringe	TA-1	\$9



FRIENDLY AND PROFESSIONAL SERVICE. CALL TODAY! 425.408.9347



MOBILE & PORTABLE DENTAL EQUIPMENT



Find The Right Portable Equipment For You

Driven to advancing dental care and patient access to care. **Call us today! 425-408-9347.**



	GO Ultra Portable	Express Portable	Transport III Portable	Lightweight Portable	AMC-15 Mobile Dental Cart	AMC-25 Mobile Dental Cart
	Our most advanced technology	World's #1 selling portable system	Military-grade performance	Only 17 pounds! Requires Portable Compressor	Affordable workhorse for general dentistry	Highly customizable, Hospital-grade dental cart
Recommended Use						
General Dentistry	✓	✓	✓	✓	✓	✓
Dental Hygiene	✓	✓	✓	✓	✓	✓
Dental Sealant Placement	✓	✓	✓	✓	✓	✓
Rotary Endodontics			✓			Optional
Motor and Handpiece Connection Type						
	Electric	Air	Electric	Air	Air	Air with Electric Option
Special Features						
Wireless Foot Control	✓					
Handpiece Illumination	✓	Optional	✓	Optional	Optional	Optional
Integrated Ultrasonic Scaler Option	AEU-350S		AEU-525S			Optional
Integrated Curing Light Option						Optional
Weight						
	32 lbs.	43.5 lbs.	56 lbs.	17 lbs.	139 lbs.	167 lbs.
	Qualifies as a carry-on luggage	Can be checked as luggage	Extra luggage fee required	Check with your airline	Usually moved by truck or van then rolled room-to-room in institutional settings	
General Pricing*						
	\$8,800-\$10,600	\$6,300	\$12,900-\$14,800	\$1,990+	\$13,500	\$16,800

Taskforce Portable Dental System

ADU-12DCE

- New look and case design
- Lightweight and easy to carry, under 33lbs!
- Two air handpiece connections
- Powerful saliva ejector vacuum
- Built-in maintenance free compressor
- Water port for Cavitron/scaler
- 110V or 220V switchable voltage
- 12.5" x 17.5" x 7.5" (31.8 x 44.5 x 19.1 cm)

The most compact pneumatic portable dental system, the Taskforce, represents the best value on the market. At 32.6 pounds the Taskforce is easily carried from location to location. This unit includes two air-driven handpiece connections (high speed /low speed), low volume suction, and air water syringe, making it the ideal system for mobile hygiene. Simple and easy to use, the Taskforce can be set up or put away in just minutes.

Taskforce **ADU-12DCE** \$4,250

Optional Accessories

- Scaler Titan Blis-Sonic Air Operated Swivel Type **ASC-06** \$1,449
- Scaler Titan Blis-Sonic 4-Hole To Swivel Adapter **ASC-06-SW** \$373
- Scaler Titan Blis-Sonic Air Operated 4-Hole **ASC-06-4H** \$1,509
- Stainless Steel Cuspidor with Solids Trap 5" (12.7mm) diameter for convenient expression during treatment **AA-26SS** \$70
- High Volume Tip Adapter **AA-12HVE** \$55 enables use of standard 11mm high volume tips
- Fiber Optic Upgrade Kits**
- 5-hole **AA-19LED-04** \$720
- 6-pin **AA-19A-06** \$720
- 4 Hole Fiber Optic Coupler **AHP-82EGEN-QD** \$295

"I have used the Aseptico Taskforce Deluxe in 25 developing countries on dental missions. The unit is lightweight, and has always functioned well with minimal power requirements. I have recommended its use to many needing a lightweight, reliable dental operating unit."

COL (Ret) Robert Meyer, DMD
Executive Director
Christian Dental Society



4-Hole to E-Type Air Motor AHP-50
Low Speed Air Motor (4-hole) adapts air turbine to accept E-type handpieces.

Call for more information.



4-Hole Restorative AHP-96MK
Air-driven
Internal irrigation
Non-optic
Max 320k rpm
Max volume 62 dB

Recommended for high speed dentistry
See page 16.



ASEPTICO PORTABLE SYSTEMS ARE MADE IN THE USA. CALL TODAY! 425.408.9347



"It's really user-friendly and quiet. Everything takes half the amount of time to do which allows me to spend more time with my patients. I've never worked with a company that cares so much about keeping me in business so that I can help take care of other people."

Erin Poicus, RDH
Cascade Hygiene

Perfect with the NEW GO Cart
Learn more about the GO cart on page 34.

GO Ultra Portable Dental Delivery System

Our Most Advanced Technology

AEU-350/AEU-350S

- Simple to operate for restorative & hygiene
- Integrated brushless motor with LED illumination
- Powerful vacuum and air flow
- Exceptionally quiet operation
- Ultra-compact and lightweight
- Meets carry-on luggage requirements
- Air filter reduces airborne particulates
- Use anywhere in the world with electric or battery power
- Customer-replaceable bacterial-viral filter
- Health Canada Approved
- Now CE-Marked (EU-MDR)

Designed for Ease-of-Use and Reliability
The GO leverages Aseptico's expertise in crafting leading portable dental units for today's challenges and opportunities. Designed for operator ease-of-use and outstanding reliability, the GO is capable of restorative dentistry and hygiene functions.

At only 32 pounds, the GO Ultra Portable dental system is so compact it qualifies as carry-on luggage for air travel.

The GO delivers powerful airflow and suction yet operates more quietly and cleanly than standard air compressor units (it even has an air filter to reduce airborne particulates).

NEW GO With			
Wireless On/Off Switch	AEU-350-W	\$9,210	
NEW GO With Scaler & Wireless On/Off Switch	AEU-350S-W	\$11,010	
GO Ultra Portable Dental Delivery System	AEU-350	\$8,800	
GO With Scaler	AEU-350S	\$10,600	

Options

24 VDC Battery	AA-350-BATCAB2	\$110
Power Cord		
Vacuum Hose Plug	AA-350-VPLUG	\$90
Snap-On Tray	AA-350-TRAY	\$110

Replacement Parts For AEU-350S

Thin Subgingival	ASC-10-PE37	\$140
Slim Universal	ASC-10-PE38	\$140
Power Universal	ASC-10-PE39	\$140
Scaler Holder	AA-350-HLD1	\$34
Scaler Wrench	ASC-10-UNIV	\$55
Scaler Ergo Sleeve	ASC-10L-ERGOSLV	\$55

Replacement Parts

Micromotor with 5' cable	AE-242SC-40	\$1,400
Water Supply Bottle	AA-350-BOT	\$46
Waste Collection Tank	330811	\$95
Vacuum Line Filter	AA-350-VLF	\$8
Air Filter	AA-350-FIL	\$43
HVE Valve and Hose Assembly	AA-350-HVE	\$180
Saliva Ejector Valve and Hose Assembly	AA-350-SAL	\$60
Foot Switch	AE-7PM	\$220
Power Cord	840107	\$23
Power Cord Pouch	410222	\$14
Hose Storage Pouch	410225	\$4

Powerful Electric Motor

5k to 200k rpm range is ideal for prophylaxis and restorative.

Faster Drying for Sealants

- Syringe Performance:
- 80 psi initial
 - Up to 50 psi continuous

More Powerful Suction!

New vacuum blower technology

NEW Ultrasonic Scaler Option Now Available

Fully integrated piezoelectric ultrasonic scaler with LED illumination

Super Compact!

22" x 13" x 9.2" fits into the airline overhead as carry on luggage

Easy and Light
Only 32 pounds!

Quieter Operation

A quieter vacuum system that turns on only when HVE and saliva ejector are removed from their holders creates a quieter experience for both operator and patient.

500 mL Water Capacity

Plenty of water capacity for restorative cases

1 Liter Waste Capacity

Easy, clean removal of waste

Impact Resistant

Durable, wheeled case with comfortable handles

Weighted & Omnidirectional Foot Pedal

Stays in place and can be activated from any angle

Wireless Foot Control
Now Available for the GO



Plug in Anywhere in the World

100-240 VAC & 24VDC input options

Can Operate on External Battery Power

Powering the GO from a 24 VDC battery is an option for those who work off-the-grid. The cable is purchased from Aseptico, and the customer supplies the 24V battery.



OUR MOST ADVANCED TECHNOLOGY!
CALL TODAY! 425.408.9347

Portable Supply Cart GO Cart

ATC-350ZT

Move equipment effortlessly with the lightweight Aseptico GO Cart. Designed for mobile dental professionals, it offers smooth all-terrain mobility, organized storage with six pouches, and built-in water bottle holders. Weighing only 16.5 lbs, it also features a durable frame that doubles as a portable seat.

GO Cart ATC-350ZT \$805

Cart size: 24.5"H x 22"W x 14.5"D (including wheels)
(51"H with fully extended handle)



"The GO Cart has completely transformed my practice. It's incredibly durable, keeps me organized, and makes moving between patients a breeze."

Sonya Dunbar, MPH, RDH
Geriatric Toothfairy | Mobile Dental Xpress



Durable Telescoping Handle
Extends Up to 51".

Perfect With The GO
Pair the GO Cart with the GO (AEU-350/s) and bring everything you need in one trip.

Extra Storage
Use the two additional side pockets for extra items.

Patented Steel Frame
Doubles as a portable seat, supporting up to 300 lbs for on-the-spot convenience.

All-Terrain Tires
Ensure smooth transportation over any terrain, making it ideal for remote and diverse settings.



Storage For Everything
Equipped with six storage pouches and customizable organizer boxes, perfect for storing handpieces, scalers, and other dental tools.



Two water bottle holders at the back can also carry sanitation supplies

Express Portable Dental System

ADU-17A and ADU-17A2

ADU-17X and ADU-17X2

- Fast and efficient for dental sealants
- Available with one or two air handpiece connections
- High volume evacuator and saliva ejector
- Built-in maintenance free compressor
- Water port for Cavitron/scaler
- High capacity 1 gallon waste container
- 110V or 230V switchable voltage
- 44-48 lbs (20-21.8 kg)
- 10.5" x 14.5" x 24" (27 x 36 x 61 cm)

Aseptico's most popular portable dental systems, the Express Air and Express II portable dental systems are fully self-contained, lightweight, reliable performers that provide everything a dentist or hygienist needs to perform effectively in the field. These powerful units feature up to two air-driven handpiece connections, high volume and saliva ejector suction, internal water reservoir, 3-way syringe, and integrated air compressor.

Express Air – features variable rheostat foot pedal. The handpiece only operates when the variable rheostat foot pedal is depressed; handpiece operates independently from suction and 3-way syringe. **(Dentist preferred)**

Express Air one handpiece connection ADU-17A \$6,100
Express Air two handpiece connections ADU-17A2 \$6,300

Express II – features on/off electric foot switch. All components are activated by the on/off foot switch. This configuration is sometimes desirable for special needs patients or other specialized applications. Please call for more information. **(Hygienist preferred)**

Express II one handpiece connection ADU-17X \$6,100
Express II two handpiece connections ADU-17X2 \$6,300

Fiber Optic Upgrade Kits		
5-hole	AA-19LED-04	\$720
6-pin	AA-19A-06	\$720
4 Hole Fiber Optic Coupler	AHP-82EGEN-QD*	\$295

Options		
Snap-On Tray	AA-17-TRAY	\$120
HEPA Filter	AA-17-HEPA	\$680

Handpieces sold separately
Recommended handpieces see pages 16-18.



"We love the ease of this piece of Aseptico equipment! It's essential for our business to run efficiently being able to utilize everything we would have in a private office setting while providing care in long-term care facilities. No sacrifice in power or suction either!"

Liz LaJoie, RDH
Director of Clinical Field Teams
Aria Care Partners

Transport III Portable Dental System

AEU-525S and AEU-525

- High torque electric motor with restorative and endodontic programs
- Fiber optic electric handpiece connection
- High volume evacuator and saliva ejector
- Quiet, intuitive, powerful
- Built-in maintenance free compressor
- Water port for Cavitron/scaler
- 115V or 230V Dual voltage
- Optional fiber optic piezo scaler
- Optional amalgam separator kit
- 57-63 lbs (25.8-28.6 kg)
- 17.9" x 10.4" x 22.1" (45.5 x 26.4 x 56.1 cm)

The Transport III is Aseptico's most capable portable dental system. Its features include noise reduction technology, an intuitive user interface, and 40k rpm brushless autoclavable micromotor. It is easy to set up, operate, and tear down. All components fit neatly into the impact resistant carrying case with inset handle and wheels.

Delivery module and user interface panel mount externally, placing controls and instruments at a comfortable level.

Dual level fluid sensors alert the user when the waste tank is 85% full for better procedure planning. The Transport III will then shut off the compressor before the tank exceeds capacity. The externally mounted waste container and optional amalgam separator means faster and cleaner waste disposal.

Transport III with Scaler	AEU-525S	\$14,800
Transport III	AEU-525	\$12,900

Options

Amalgam Separator Kit	AA-525AMAL	\$730
Micromotor with 6ft cable	AE-240-40	\$1,750
Foot Switch with 10ft cable	AE-7PM-10	\$240

Replacement Parts


Micromotor with 5ft cable	AE-240SC-40	\$1,750
Waste Tank Assembly	330603	\$660
Foot Switch	AE-7PM	\$220
Power Cord 110V	840049	\$44
Power Cord 220V	840007	\$14
Water Reservoir (white)	730471	\$13

Replacement Scaler Tips For AEU-525S Only

Slim Universal	ASC-10-PE38	\$140
Thin Subgingival	ASC-10-PE37	\$140
Power Universal	ASC-10-PE39	\$140



*Handpieces sold separately
Recommended handpieces see pages 16-18.*

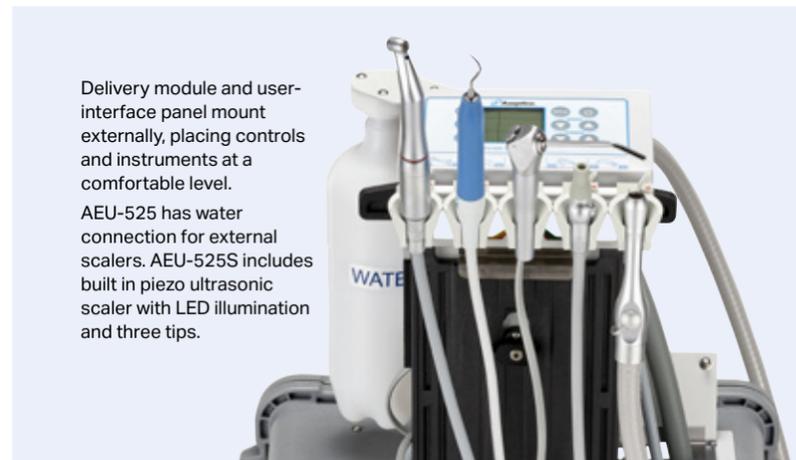

Delivery module and user interface panel mount externally, placing controls and instruments at a comfortable level. AEU-525 has water connection for external scalers. AEU-525S includes built in piezo ultrasonic scaler with LED illumination and three tips.

Lightweight Portable Delivery Systems

ADU-04ST

- Compact design and case under 15lbs (6.8kg)
- Customizable handpiece (high/slow speed) and vacuum options
- 3-way air/water syringe
- 25 ft. air line to connect to any air source
- Instrument tray
- Padded carrying case
- Carrying case: 25" x 11.5" x 11.5" (63.5 x 29 x 29cm)
- Requires air compressor for operation
- Quick disconnect options, .25 QD-Standard or 3/8 Schrader (optional)

Restorative System w/ 2 HPC Connections	ADU-04ST	\$1,990
---	-----------------	---------

Optional Upgrades

Scaler Titan Blis-Sonic Air Operated Swivel Type	ASC-06	\$1,449
Scaler Titan Blis-Sonic 4-Hole To Swivel Adapter	ASC-06-SW	\$373
Scaler Titan Blis-Sonic Air Operated 4-Hole	ASC-06-4H	\$1,509
High Volume Suction	AA-21ST	\$800
Saliva Ejector	AA-25ST	\$730
5-hole Fiber Optics	AA-19LED-04	\$720
6-pin Fiber Optics	AA-19A-06	\$720
4 Hole Fiber Optic Coupler	AHP-82EGEN-QD	\$295

Product shown with high volume suction and saliva ejector suction options added



Airpac III Portable Compressor

AA-100-110

- 1.0 horsepower, 100 psi oil-less motor
- Air delivery 3.1 CFM at 40 PSI
- Built-in wheels and pull handle
- Available in 110V only
- Not available in Canada
- Weight: 46 lbs. (20.87 kg)
- Size: 17.93"W x 22.06"L x 10.43"H (45.5 cm W x 56.0 cm L x 26.5 cm H)

The Airpac III Portable Compressor is built by Aseptico to stand-up to heavy use and travel, the compressor is secured in its own super duty carrying case with built-in wheels and pull handle. The Airpac III includes an air storage tank and automatic pressure switch control.

110V Compressor	AA-100-110	\$2,350
-----------------	-------------------	---------




ASEPTICO PORTABLE SYSTEMS ARE MADE IN THE USA. CALL TODAY! 425.408.9347



AMC-25 Mobile Dental Cart

The AMC-25 Mobile Cart is a superior solution for providing general dental care in nursing homes, hospitals, or anywhere a mobile dental unit is needed. Available optional accessories allow for customizing this dental cart to meet your needs, including the ability to add electric handpieces, handpiece illumination, and much more!

- No plumbing required!
- 2 automatic high speed air handpiece connections (4-hole; non-fiber optic)
- 1 automatic low speed handpiece connection (4-hole; non-fiber optic)
- High and low volume suction with automatic overflow protection and solids trap
- 3-way air/water syringe
- Quick connect for ultrasonic scaler
- 2 hospital-grade outlets for accessories
- Wet/dry foot control
- Air Auxiliary Quick Disconnect
- Dual voltage 110V or 220V
- Operating pressure 80 psi
- High volume vacuum 7.6 SCFM @ 1.3 in Hg, 5.6 SCFM @ 4 in Hg
- Low volume vacuum 2.2 SCFM @ 1.5 in H
- Noise Level 65 dBA or less @ 3'4" (1 meter)
- Weight 167 lb (75.8 kg)
- Dimensions 23.5" W x 30" L x 36.5" H (56.7 x 76.2 x 92.7 cm)
- Optional: integrated electric motor, fiber-optic illumination for handpieces, ultrasonic scaler, curing light, amalgam separator, external vacuum/air/water hospital connections, auxiliary arm with accessory tray

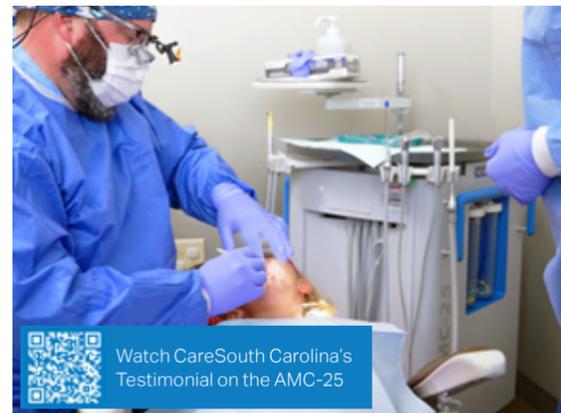

AMC-25 Mobile Dental Cart **AMC-25** \$16,800

Options

Integrated Electric Motor	330674	\$1,990
Valo Curing Light	ALU-25CLK-VALO	\$1,795
2 Fiber Optic HP Connections		
5-hole (ISO-B)	AA-19LED-25CK5	\$1,400
6-pin (ISO-C)	AA-19A-25CK	\$1,350
Piezo Scaler	ASC-10-25CKL	\$2,050
Amalgam Separator	AA-20AMAL-CK	\$920
Auxiliary Arm with Accessory Tray	AA-20A	\$520
Hospital Kit	AA-25HK	\$1,750
(external air, water and vac connections)		

Replacement Parts

Micromotor	AE-240SC-40	\$1,750
Wet/Dry Foot Control	AA-43W	\$110
Air/Water Syringe	TA-90D	\$85
1L Water Reservoir – Clear	730631-01	\$70
1L Water Bottle – Opaque (each)	730471	\$13
Strainer	730745	\$34
Amalgam Separator Filter	730595	\$160
Foot Control Caddy	AA-25CDY	\$45

Unparalleled Features & Customization

The AMC-25 Mobile Cart is not only the most advanced dental cart on the market, but it also offers a range of customizability to suit your specific needs.



Auxiliary Arm with Accessory Tray*

360 degree maneuverability to be able to operate from any angle.

Industry Leading VALO Curing Light*

Easy connect/disconnect with its own individual holder.

Scaler Connection*

Optional ultrasonic scaler connection can be hardwired in.

Multiple Fiber Optic Connections*

Optional connections for your fiber optic handpieces.

Integrated Electric Motor*

Read more about functionality on the AEU-5000 (pg. 9).



Completely Mobile for Easy Transport

Easy to wheel into trucks or vans, down hallways, in elevators, ER's and operating rooms.



Wet/Dry Foot Pedal

Stays in place and can be activated from any angle.

Amalgam Separator

Filters out particulate matter from waste water.

Foot Control Caddy

Easily store the foot control and other accessories on your AMC-15/25



1L Water Reservoir & Bottle

Large capacity for minimal water recycling between patients.

Great Backup Unit for Offices & Hospitals

Impact Resistant

Durable, wheeled case with comfortable handles.

Hospitable Grade Outlets

For any additional accessories needed.



*Does not come standard with unit



ASEPTICO MOBILE DENTAL SYSTEMS ARE MADE IN THE USA. CALL TODAY! 425.408.9347

AMC-15 Mobile Dental Cart

Designed for Easy and Affordable Mobility, the AMC-15 Mobile Dental Cart is equipped for general dentistry with two automatic high speed air connections, one low speed handpiece connection, internal vacuum and air, water supply, large waste tank with automatic overflow protection, and amalgam separator.

- No plumbing required!
- 2 automatic high speed air handpiece connections (4-hole; non-fiber optic)
- 1 automatic low speed handpiece connection (4-hole; non-fiber optic)
- High & low volume suction with automatic overflow protection and solids trap
- 3-way air/water syringe
- Quick connect for ultrasonic scaler
- Amalgam Separator
- Wet/dry foot control
- Air Auxiliary Quick Disconnect
- Available in 110V or 220V
- Operating pressure 80 psi
- High volume vacuum 6.0 SCFM @ 0 in Hg, 3.6 SCFM @ 4 in Hg
- Low volume vacuum 1.3 SCFM @ 1.5 in Hg
- Noise Level 63 dBA @ 3'4" (1 meter)
- Weight 139 lb (63 kg)
- Dimensions 23.5"W x 30" L x 36.5" H (56.7 x 76.2 x 92.7 cm)
- Optional: auxiliary arm with accessory tray

AMC-15 Mobile Dental Cart **AMC-15** \$13,500

Options

Auxiliary Arm with Accessory Tray	AA-20A	\$520
4 Hole Fiber Optic Coupler	AHP-82EGEN-QD*	\$295

Recommended Handpieces

Air-Driven see pages 16-17

Replacement Parts

3-Way Air/Water Syringe with Tip	TA-90D	\$85
Amalgam Separator Filter	730595	\$160
Clear 1-liter Water Bottles w/Level Indicators	730631-01	\$70
Two 1-liter Water Bottles	730471	\$13
Wet/Dry Foot Control	AA-43W	\$110
Power Cord 115V	840049	\$44
Foot Control Caddy	AA-25CDY	\$45



Compricart II Delivery System

ADU-20B

- One high speed and one low speed handpiece connection
- High volume evacuator and saliva ejector
- Large stainless steel instrument tray (13.5" x 9.75")
- All day water system with two liter self-contained water bottle

Roll the Compricart II Dental Cart from room to room with ease! Great for hospitals and expanding dental offices, the simple controls of the Compricart II make operation stress-free. A compact design, ample work surface, stainless steel instrument tray, integrated water supply, and venturi suction (high volume and saliva ejector) round out the features of this versatile platform. Just connect to an air source – no central vacuum required! Including the optional AWS-1C allows an external scaler to use the cart's water supply.

Compricart II	ADU-20B	\$3,750
Options		
Water Connection for External Scalers	AWS-1C	\$85
Fiber Optic Upgrade Kit 5-hole	AA-19LED-04	\$720
Fiber Optic Upgrade Kit 6-pin	AA-19A-06	\$720
Replacement Parts		
Complete HVE Assembly	330461	\$540
Disc Foot Control w/ Water Toggle	AA-43W	\$110
Air/Water Syringe	TA-90D	\$85
2L Water Bottle	730665	\$16



AEU-5000 Add-on Option Available
See details on page 9

"We have been using our AMC-15 since November 2023 approximately once per month for dental care in people with severe developmental disabilities. It is the perfect unit for our needs and has functioned flawlessly."

Gregory A. Costello, MD

Chief Medical Officer
Hope Christian Services



FRIENDLY AND PROFESSIONAL SERVICE. CALL TODAY! 425.408.9347



NEW Sileather Fabric

Discover Sileather™ — eco-friendly and sustainable silicone fabric designed for durability and easy care. Naturally stain-resistant, UV-stable, and tough against wear, Sileather keeps our stools and chairs reliably clean and comfortable in demanding settings.



Chairs and Stools Available in Gray and Blue



All seating includes a padded carrying case



AseptiStool

ADC-08

- Set-up and store in seconds
- Includes padded case – no loose parts to misplace!
- 300 lb weight capacity
- Adjustable seat height: 19" - 22" (48.3-55.9 cm)
- 15 lbs (6.75 kg)
- Padded carrying case: 18" x 22" x 8" (48 x 53 x 23 cm)

Durable, adjustable and foldable, the AseptiStool Portable Dental Stool is the perfect seating solution for any clinician. Its one-piece design sets up or stores in seconds, and it swivels on four casters for easy mobility during treatment. It is height adjustable, lightweight, and folds up for simple storage in the included padded travel case.

AseptiStool	ADC-08	\$1,020
Replacement Padded Carrying Case	410093	\$190



Hydraulic AseptiStool

ADC-10

- Easy set-up and storage
- Includes padded case
- 300 lb weight capacity
- Hydraulic adjustable height: 18" - 25.5" (45.7 - 64.8 cm)
- 19 lbs (8.6 kg)
- Padded carrying case: 26" x 26" x 11" (56 x 56 x 28 cm)

With hydraulic adjustable height, the Hydraulic AseptiStool gives the mobile dentist the comfort and flexibility of doing dentistry in the office. It is easy to set up, tear down and transport in its included padded case. Even with its sturdy construction and high weight capacity, the stool is lightweight and carrying it into any mobile dental scenario is simple.

Hydraulic AseptiStool	ADC-10	\$1,350
Replacement Padded Carrying Case	410213	\$250



AseptiChair

ADC-01

- Sturdy, high-strength steel construction
- Lightweight, only 30 lbs (13.5 kg) in case!
- 300 lb weight capacity
- Backrest adjusts from fully upright to supine
- Price includes padded case, arm slings, and headrest cushion
- Adjustable seat height: 11" - 22" (11 - 55.9 cm)
- Padded carrying case: 23" x 34" x 8" (58 x 86 x 20.3 cm)

The AseptiChair's sturdy rolled steel construction is lightweight yet strong, making it perfect for portable dentistry. It is adjustable to five positions, with seat height ranging from 11 to 22 inches. The toe board folds down, allowing this portable dental chair to be used for examination.

AseptiChair	ADC-01	\$2,200
-------------	---------------	---------

Replacement Parts

Replacement Padded Carrying Case	410193	\$290
Replacement Leg Adjustment Pin	460445-08	\$17
Endcap Heavy Duty 1"	850026	\$3
Armrest Support Pin	460444-08	\$9



Hydraulic AseptiChair

ADC-02

- Hands-free height and back adjustment
- Adjustable with patient in chair
- 300 lb weight capacity
- No external electricity or air required
- Handle and wheels for easy mobility
- 69 lbs (31.3 kg)
- Hydraulic adjustable seat height: 11" - 33" (27.9 - 83.8 cm)
- Folded dimensions: 33" x 22" x 18" (83.8 x 55.9 x 45.7 cm)

The Hydraulic AseptiChair provides hands-free operation and is the only portable chair available with adjustability functions comparable to a standard dental chair. All chair lift functions are foot actuated hydraulic for hands-free operation; no external power is required. The AseptiChair provides comfortable patient seating for patients up to 300 lbs. Easy move wheels.

Hydraulic AseptiChair	ADC-02	\$7,550
Hydraulic Assembly Only	330629	\$2,100

Replacement Parts

Replacement Dust Cover	410210	\$420
------------------------	---------------	-------

"Having quality products like Aseptico patient/operator chairs gives our providers, patients, and facilities confidence in the care that a mobile alternative practice can provide"

Sara Luetke, RDH
Sound Dental Care



ALL ASEPTICO SEATING IS MADE IN THE USA. CALL TODAY! 425.408.9347

NOMAD™ PRO 2 Defining Portable Dental X-ray Technology

ARU-07KP2

- Best-in-class portable digital X-ray system
- Lightweight 5.5lbs (2.5kg)
- Simple touch controls on color LCD screen
- Removable and rechargeable battery packs
- One battery charge allows for more than 100 exposures
- Heavy Duty Protective Transport Case Included



The NOMAD Pro 2, Portable Dental X-ray System captures clear radiographic images with ease – anytime, anywhere. Its lightweight, ergonomic, handheld design makes it ideal for mobile dental units, or practices with multiple operatories. Best-in-class features include advanced shielding to ensure operator and patient safety while delivering precise imaging results.



Offering cordless, handheld operation, the NOMAD Pro 2 portable x-ray system is powered by rechargeable battery packs and comes contained in an easily transportable, rugged case with pull-out handle and wheels. Long-lasting batteries deliver exceptional performance, allowing you to complete hundreds of exposures on a single charge. Intuitive user interface and quick setup save valuable time, enabling you to focus on what matters most—providing patient care.

Nomad Pro 2 Portable X-Ray System with Carrying Case **ARU-07KP2** Call For Best Price
Includes Premium Upgrade Kit

Replacement Parts

Premium Upgrade Kit
Includes Two Replacement Battery Handsets
One Replacement Charging Cradle
Power Cord **ARU-10134502** \$953

Basic Upgrade Kit
Includes Replacement Battery Handset
One Replacement Charging Cradle
Power Cord **ARU-10135062** \$589

Replacement Battery Handset **ARU-10131887** \$519
Replacement Charging Cradle w/Power Cord **ARU-10132534** \$270
Not Available In Canada



Portable LED Exam Light

ALU-40LED

- Adjustable LED with seven beam diameters
- Durable construction
- Includes padded carrying case
- Optional Amber UV filter for working with composites
- 12 lbs (5.4 kg)
- Padded carrying case: 45" x 9" x 5.5" (114.5 x 23 x 14 cm)

The Portable LED Exam Light provides high-intensity LED light with seven beam diameter settings. It features a compact lamp head with adjustable intensity, flexible gooseneck, telescoping tripod stand and padded carrying case. When a standard electric outlet is not available, an optional auxiliary DC input can be purchased that allows 12V and 24V battery operation in the field.

Portable LED Exam Light	ALU-40LED	\$1,630
Optional Amber UV Filter	461974	\$4
Replacement Padded Carrying Case	410084	\$180
DC Cable with Battery Clips	AA-40DC-BC	\$170
DC Cable with Car Lighter Adapter	AA-40DC-LP	\$200


Portable Tray Stand

ATC-03CF

The Portable Tray Stand is a sturdy and necessary addition to any mobile operatory. Standard size stainless instrument tray removes easily for autoclaving.

It weighs less than 9 pounds and folds into included case when work is done. Made in USA.

Portable Tray Stand **ATC-03CF** \$990

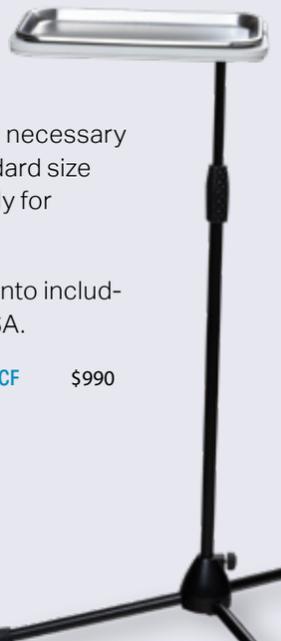

CALL TODAY! 425.408.9347

FOR FRIENDLY AND PROFESSIONAL CUSTOMER SERVICE.

Heavy Duty Brushless Lab Motor Torque Pro Delivers Power

AEU-16SS-120

- For heavy acrylic use
- Brushless 50k rpm micromotor
- Consistent torque even at low speeds
- Includes Variable Speed Foot Control



The Heavy Duty Brushless Lab Motor delivers power, performance, durability, and portability. The brushless 50k rpm micromotor can produce up to 8.3 Ncm of torque and features easy to read digital display, preset memory storage, and hand/foot change button.

Heavy Duty Dental Lab Motor **AEU-16SS-120** \$1,200

Replacement Parts:

Motor/Cable Assy	AE-H180	\$540
Variable Speed Footswitch	AE-SP	\$80

Light Duty Lab Motor Powerful & Lightweight

AEU-10SS

- For light to medium acrylic use and low speed prophy
- Versatile E-type straight 1:1 handpiece included
- 20k rpm micromotor
- Available in 110V or 220V
- Includes AE-7P On/Off Foot Switch

The Light Duty Lab Motor is one of the finest motor and handpiece combinations available to the professional. It is ideal for chairside adjustments, prophylaxis, or any light to medium use.

Light Duty Dental Lab Motor **AEU-10SS** \$1,150

Light Duty Dental Lab Motor **AEU-10SS-220** \$1,200

Options

Optional Variable Speed Foot Control **AE-14A** \$360

Replacement Parts:

1:1 Handpiece	AHP-101	\$380
Micromotor	AM-103B	\$230
On/Off Foot Switch	AE-7P	\$220



Inludes AHP-101 Handpiece

Battery Powered Portable Lab Motor Ideal for Mobile Hygiene

AEU-03M1B

- For light duty use and low speed prophy
- Straight 1:1 e-type handpiece included
- Dual voltage 110V-240V

The Battery Powered Portable Lab Motor is a smart solution for lightweight lab benchwork and performing prophylaxis in the field! This small portable unit delivers 22k rpm and can perform 8 hours on a single battery charge.

Portable Lab Motor with Handpiece **AEU-03M1B** \$850

Portable Lab Motor **AEU-03M1B-NHP** \$600

Replacement Parts

1:1 Handpiece	AHP-101	\$380
Micromotor	AE-5M1-25	\$180
Foot Switch	AE-7M1	\$23



Portable Sterilizer

ASU-12-110 and ASU-12-230

- Compact and portable
- 18-minute sterilization cycle (24-minute total cycle time)
- Only 6 minutes to reach sterilization temp
- Temperature: 250°F (121°C)
- 11.5 lbs (5.2 kg)
- 9 liter capacity



The Portable Dental Autoclave is designed to sterilize solid, unwrapped instruments in a 121°C saturated steam process during 18-minute sterilization cycle. Use distilled or deionized water only.

110V Sterilizer **ASU-12-110** \$1,450

230V Sterilizer **ASU-12-230** \$1,550

Portable Saliva Ejector

AA-20

- Great for mobile dentistry or as a backup for sedation dentistry
- Powered by rechargeable battery pack, wall outlet or 12 VDC vehicle source
- Up to 60 minutes of continuous use on a single charge
- 800cc waste canister with float shut-off
- 6.3 lbs (2.9 kg)
- 9" x 7" x 8" (22.8 x 17.8 x 20 cm)



Great for mobile dentistry or as a backup for sedation dentistry, the Portable Saliva Ejector can be powered by a wall outlet, 12V DC source, or rechargeable battery pack. It features an adjustable 80 - 550 mm Hg vacuum range and 800cc waste canister with float shut-off. Includes padded carrying case.

Portable Saliva Ejector **AA-20** \$650

Replacement 800 ml Canister **AA-20CAN** \$40

Bacteria Filters **AA-20BF-2** \$12

Vacuum Pump Tube Assembly **330668** \$49

In addition to providing easy portable suction, the Portable Saliva Ejector fulfills the legal requirement for a back-up suction system in the event of power or suction failure.

Shop online at aseptico.com



CALL TODAY! 425.408.9347

FOR FRIENDLY AND PROFESSIONAL CUSTOMER SERVICE.

Customer Service 425.408.9988

Monday-Friday
6:00 a.m.- 3:30 p.m. PST
Customer service representatives
are standing by to help.
repairs@aseptico.com

Quality System

Aseptico's Quality Management System is certified to ISO 13485 as well as the Medical Devices Single Audit Program Recognized By The U.S. and Canada.

Return and Repair Policies

- Products that have been in clinical use cannot be returned for credit.
- Service and repair charges are billed on credit card. Customer is expected to accept charges up to \$750 without advance notification. Minimum service/repair charge is \$160.
- All returns and repairs require prior authorization. Contact Aseptico Customer Service for a Return Material Authorization number (RMA#).
- The RMA# must appear clearly on the exterior of the package being returned or it may be rejected and returned to the sender.
- Packages received without prior authorization are not the responsibility of Aseptico. Aseptico reserves the right to dispose of such packages and their contents 30 days after receipt with no obligation to the sender.
- The sender is responsible for proper packaging and for all costs for return shipment to Aseptico.
- New or unused products may be returned for full credit (less shipping charges) under a Return Material Authorization number within 30 days of date of original invoice.

Warranty

Aseptico warrants its products against defects in material or workmanship for a period of two (2) years, from date of original invoice. This warranty is not transferrable. Some handpieces are warranted for one year under the same conditions. Other handpieces and expendable components, such as air turbines and light bulbs, are covered by shorter warranty periods, or have no warranty. Aseptico's sole obligation under product warranty is (at its sole option and discretion) to repair or replace any defective component or product in part or whole. Aseptico shall be the sole arbiter of such action. In the event of alleged defect under warranty, the purchaser is to notify Aseptico's Customer Service Department promptly. Customer Service will provide instructions, usually directing that the product be returned for service. Shipment to Aseptico and the cost thereof is always the responsibility of the purchaser.

Accidental misuse, inappropriate installation, or failure to perform directed maintenance voids the warranty.

Aseptico does not assume, under this warranty, any risks or liabilities arising from the clinical use of its products, whether or not such use involves coincidental utilization of products manufactured by others. Deliberately defacing, modifying, or removing the serial number voids the warranty.

Sales 425.408.9347

Monday-Friday
6:00 a.m.- 4:30 p.m. PST
Sales representatives
are ready to assist you.
orders@aseptico.com

Ordering Policies

- Aseptico accepts Visa, MasterCard, Discover, and American Express.
- Open account is available only if credit is approved in advance.
- Freight costs are added to all invoices.
- All orders are shipped FOB origin.
- Delinquent accounts are subject to 1.5% monthly interest, attorney's fees, collection and court costs when needed to satisfy debt.
- Prices subject to change without notice.
- Please visit: aseptico.com/warranty-policies for additional warranty and repair policy information.

Sales Tax

Aseptico uses a tax service (Avalara) to calculate and apply state sales tax to invoices for all states that Aseptico is registered to collect/remit sales tax. We reserve the right to register in additional states at any time and begin collecting/remitting sales tax.

The amount of tax charged on your order will depend upon many factors including type of item purchased and destination of the shipment. Factors can change between the time you place an order and the time of invoicing, which could affect the calculation of sales taxes. The amount appearing on your order as sales tax may differ from the sales taxes ultimately charged. Customer is responsible for any applicable state sales or use taxes not charged by Aseptico.

Sales Tax Exemptions

Unless Aseptico is in possession of a validated exemption certificate prior to shipment, the customer must pay the entire invoice, including sales tax. Certificates can be emailed to AR@aseptico.com. You can upload your Tax Exempt certificate using our Avalara Cert Express Tool: Aseptico.com/tax-exempt. Please inform your sales rep when you do so, and they will get it applied to your account for future orders placed through your sales rep.

International Customers

Most Aseptico products are CE-marked and Health Canada licensed. However, customers should verify the status of any specific item(s) with their sales representative or local distributor prior to finalizing the purchase. Please contact your local Aseptico distributor for pricing and availability.

Prices shown in the catalog are in US dollars. If you are unsure of the Aseptico distributor in your country, for specific ordering and payment policies, or for general product inquiries, please contact our International Sales Department at int@aseptico.com.

Aseptico, Inc.
8333 216th St. S.E.
Woodinville, WA 98072

The Aseptico Global Dental Catalog is published annually by Aseptico Inc., 8333 216th Street SE, Woodinville, WA 98072. Aseptico is a global leader in dental equipment design and engineering. The company has endeavored to represent all current product information, specifications and pricing accurately. Product information may change. If you would like to communicate with Aseptico, the company welcomes your interest and invites you to contact us at 425-408-9347 or via aseptico.com. ©Aseptico Inc., 2026 all rights reserved. "Aseptico" is a registered trademark of Aseptico Inc.



Aseptic, Inc.
8333 216th Street SE,
Woodinville, WA 98072

20% OFF 7000 Series Motor Package

AEU-7000L-70V-W

Only \$5,112 Save \$1,278

AEU-7000L Motor with LED illumination, Mont Blanc fiber optic handpiece, and the new Wireless Variable Foot Control

Call 425-408-9347 To Learn More



Portable Dentistry GO Package

Purchase a GO with wireless foot pedal (AEU-350-W or AEU-350S-W) and get a FREE GO Cart (ATC-350ZT) and Highspeed Handpiece (AHP-72MK-FO)

Call 425-408-9347 To Learn More



DENTAL MOTORS **WE POWER DENTISTRY**

MOBILE and PORTABLE DENTISTRY **MOVING TO MEET NEED**

DENTAL EQUIPMENT OEM **EXCELLENCE IN ENGINEERING**

GLOBAL DENTAL CATALOG **PROVEN PRODUCTS**